

尺寸:670*150mm 材质250克铜板纸 正面光膜 彩色印刷 (光标外为裁切部分)

DeenSmart
Smart & Safe


USER MANUAL

Wireless Massage Pads



P30

Wireless Massage Pads Product Descriptions

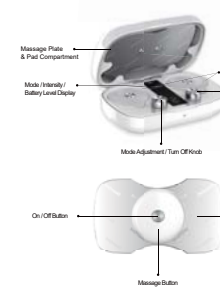

Package Content

- 1 x P30 Wireless Massage Pads (120*79*39 mm)
- 2 x Massage Button
- 4 x Massage Pad (100*50 mm)
- 1 x Massage Pad Storage Plate
- 1 x USB Cable
- 1 x User Manual

3. The magnetic gel pads on film is consumable and can be used 20-40 times. Its service life varies according to washing frequency, skin condition, and storage state. When viscosity drops, please consider purchasing a new magnetic gel pad. Please contact the retailer to purchase the magnetic gel pad of the same model and specification.

Troubleshooting

Fault Phenomenon	Cause Analysis	Resolution Method
Message Pads cannot turn on.	Battery power is exhausted.	Please connect the adapter output standard 5V1A to charge.
No pulse output (no stimulation or weak stimulation when the magnetic gel pad is attached to the skin).	The magnetic gel pad is low power. The magnetic gel pad is damaged or contaminated, not sticky, and has poor contact with the skin.	Please connect the adapter output standard 5V1A to charge. Replace with new magnetic gel pad.
There is less strong feeling on the skin.	Used too long. The intensity of use is too great. The magnetic gel pad does not stick to the skin firmly. Skin allergy to magnetic gel pad.	The appropriate time is 15-30 minutes. Decrease the intensity properly. Please stick the magnetic gel pad firmly stick to the skin. If you have a history of skin allergy, please stop the use only allergy, please attention the message time of the same skin area.
The remote control is not working.	Please check if the magnetic gel pad is dirty. The message button is beyond the control range of the remote control. When the message button is placed in the desk, the position is not aligned. The charging pin does not press well.	Wipe the skin with a damp towel. The appropriate time is 15-30 minutes. Decrease the intensity properly. Please stick the magnetic gel pad firmly stick to the skin. If you have a history of skin allergy, please stop the use only allergy, please attention the message time of the same skin area. Clean the use to replace the magnetic gel pad film. Please keep the remote distance of the message button and the control box within 5-10 meters. Please make the message button not necessarily to align the charging pin at the desk. When the message button is placed in the desk, the position is not aligned. The charging pin does not press well.
Cannot charge.	The base cannot be charged, the charging cable is not good. The message button is damaged, the port is not charged properly or in poor contact.	Replace the USB charging cable. Check if the port is plugged in properly or replace the power cable.

Note: When the service life of the message pads ends, please discard the message pads built in the lithium battery in an environmentally harmless treatment according to the local laws and regulations.

Mode Descriptions

Mode 1	Body Massage Used for daily relaxation, the massage helps to relieve fatigue and refresh the body.	
Mode 2	Neck Massage Used for neck health, the massage improves blood circulation of the neck and relieves stiffness and discomfort.	
Mode 3	Shoulder Massage Used for shoulder relaxation, the massage helps to relieve shoulder stiffness and discomfort.	
Mode 4	Wrist Massage Used for wrist health, the massage improves blood circulation of the wrist and relieves stiffness and discomfort.	
Mode 5	Low Massage Used for low back health, the massage improves blood circulation of the low back and relieves stiffness and discomfort.	
Mode 6	Joint Massage Used for joint health, the massage improves blood circulation of the joint and relieves stiffness and discomfort.	

Important Safety Instructions

- The message pads is not a medical device and cannot replace medicine for pain treatment.
- People with the following medical conditions implanted in the body, the use of message pads is strictly prohibited, such as cardiac pacemakers, artificial lungs, etc.
- If you are receiving medical treatment from your doctor, please consult with your doctor before using this message pads.
- The skin areas are with wound, healing, or inflammation on the surface of the skin, please stop using this message pads during this time.
- People with sensitive skin may be allergic during use. There is a small release under or around the message pads, and it usually disappears in about 2 hours. If the body has persistent abnormality or serious discomfort, please stop using it at once and consult your doctor immediately.
- People with sensory impairment in the skin, children and people who cannot fully express their wishes are strictly prohibited from using this message pads.
- Do not use the message pads on the breast (mammary gland), on the head, or the heart, and avoid any sensitive areas of the body.
- Do not use the message pads on people with metal implants, including x-ray equipment, metal, in case of electrical shock.
- Do not use the message pads on wet patients.
- Do not use the message pads during pregnancy.
- Do not use the message pads on the genital area.
- Please keep this message pads out of the child's reach to avoid children swallowing or damaging message pads.
- When the message button is out of the box but not connected to the message pad, please keep it in a safe place to avoid the finger touching the charging pins on the message button to produce a micro current. This feeling is harmless, but it will cause discomfort. And do not remove the message pads while message is avoid touching discomfort.
- Do not use this message pads while bathing, otherwise, it may cause electric shock or short circuit.
- Do not use the message pads in bathroom, kitchen, and electromagnetic radiation places.
- Do not use the message pads while sleeping, and do not stick this message pads the same area of skin for a long time without taking off to cause skin irritation.
- The same person be charged, the port is not charged properly or in poor contact.
- Do not repair or disassemble this message pads by yourself.
- Do not use in specific places such as hospital, large equipment, swimming pools, bathrooms, etc. in case of danger.
- When the product is discarded, please do not discard the lithium-ion battery in the product at will, please dispose of it in accordance with the relevant local laws and regulations.

About This Product

Thank you for purchasing our Wireless Massage Pads. This is a stylish and innovative massage solution, with multiple modes and intensities, that offers a convenient portable and wireless operation. It has a built-in high-quality long-life standby lithium battery, using Electrical Muscle Stimulation (EMS) and Transcutaneous Electrical Nerve Stimulation (TENS) technology. The integrated design of the charging case has multiple functions, including remote control for the selection of modes and intensity, along with built-in storage compartments for the message pads and wireless charging, with the convenience of portability in a sleek and modern case.

The Wireless Massage Pads include various modes for personal massaging techniques, they can promote local blood flow, and effectively relieve fatigue, stiffness, numbness, aches, and other discomforts of the body's neck, shoulder, back, limbs, and joints.

Our Wireless Massage Pads can meet different needs of different people, including daily health care, working fatigue care, and muscle relaxation after exercise, which is worth of choice due to its wide proposition. Use our Wireless Massage Pads anytime, anywhere, whenever you want them, and use them freely to take care of your body and health.

Product Features

- The message pads, multi-function integrated design, utilizing wireless remote control technology, high definition display, knob control operation, simple, easy, and convenient to use.
- The message pads control case has built-in storage and charging functions, each Pad has independent and automatic charging, with a long battery life that will secure constant operation.
- Equipped with dual buttons, massage different areas at the same time. Storage of pads and pads, a conveniently way.
- Professional massage techniques, using 6 different modes and a 10-level of intensity adjustments.

Operation Instructions

- Open the box, take out the message pads, place them horizontally, and align the pins to connecting the pins on the message button, it will just fit through a magnetic connection, and then use off the transparent protective film on the message pads.



Product Technical Specifications

Product Name: Wireless Massage Pads
Model: P30
Voltage: DC5V1A
Built-in Lithium Battery
Control Box: 3.7V1000mAh
Message Button: 3.7V1000mAh

Normal Working Conditions:
Ambient Temperature: 5°C-40°C
Relative Humidity: 10%-90%
Atmospheric Pressure: 950Pa-1050Pa

Electrical Performance:
Pulse Frequency: 0Hz-120Hz
Pulse width: 20us-600us
Rated Load Impedance: 500Ω, Tolerance ±10%
Message Working Time: 15 minutes, Tolerance ±10%
Message Button Life: 5 years
Battery Life: 100 times charging

DEEN SMART WARRANTY

- Guarantee: If the product is under normal use and maintenance (follow the manual instructions), due to material, workmanship or manufacturing problems resulting in performance failure, during the warranty period, we will provide you with a replacement service and will not charge you for the warranty service. If the retailer has an additional promise to you, you should ask the retailer to honor it, and Deen Smart is not responsible for this.

- Warranty Period: One year from the purchase date of the product (subject to the invoice date of the purchase of the product).

- Warranty Claim Process: Customer must fill in the form (provided by Deen Smart Service Team) with a detailed description of the issue and sufficient proof of purchase (the order number from our authorized retailers, pictures, and videos of the defective product) (the use of any other product inquiry, all unsubstantiated claims will be rejected).

- Note: It may be necessary to return a defective product for quality inspection before Deen Smart releases a product replacement.

- Not Covered Under Warranty
 - Purchased from unauthorized retailers or third parties.
 - Products without sufficient proof of purchase.
 - Items that have expired their warranty period.
 - Non-quality related issues.

- Use a clean slightly wet towel (damp but not dripping water), and wipe off the sweat on the skin surface where the pads will be placed.



- Stick the message pads firmly on the skin, and long press the switch button in the middle of the message button 5-10 seconds when it is turned on. The buzzer will have prompt tone sounds and the indicator will light up. After turning it on, if the message button detects that the message pad is not making direct contact, the message button will turn off after five times, entering into the power-saving mode.



- The message pads can be affixed to different parts of the body respectively, and they can be operated and controlled simultaneously by the wireless remote control to achieve a synchronized massage.



- When opening the control box it will be in standby mode, the operation panel lights up in default mode, and the power indicator is turned on. If the control box light button, it indicates low power mode and it will need to be charged. The mode and intensity knobs on the operation panel are adjusted with a buzzer sound to freely select the message mode (1-6) and message intensity (01-10). Each mode has 15 minutes of operation, and after that, it will automatically be turned off.



- A. The mode knob, increases the intensity clockwise, decreases the mode counter-clockwise, increases or decreases the mode step by step, there are 6 modes in total. When adjusting to mode 6 or mode 1, the user is adjusting the knob will be with a buzzer sound, generally indicated each time, the buzzer prompt once with the rotate knob. The mode knob, press up and down, also has the function of shutdown to stop dual pads massaging simultaneously for convenience.

FCC Warning

- Free products.
- Repairs through 3rd parties.
- Damage from misuse of products (including, but not limited to: falls, operating devices improperly).
- Failure or damage caused by force required or unexpected events.
- Product repairs, rust and damage caused by contact or exposure to improper temperature, solvents, acids and alkalis, prolonged immersion in water or foreign matter intrusion, etc.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If the equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

Product Information

Please visit our website and social networks for more information about our Wireless Massage Pads:

- www.deen-smart.com
- www.facebook.com/deensmart
- www.instagram.com/deensmart

Contact information

You can also contact us via email with any other product inquiry

support@deen-smart.com