

# ROBINT SMART BOX

## USER MANUAL



ACCOUNT QR CODE



Beijing Robint Technology Co., Ltd.

Website: [www.irobint.com](http://www.irobint.com)

Address: Room 2006, Block A, Jinqiu International Building,  
Zhichun Road, Haidian District, Beijing.

Customer service phone number: 010-

BEIJING ROBINT TECHNOLOGY CO., LTD.


## ABOUT PRODUCT SPECIFICATION

The copyright, trademark and trade name of the instruction manual are owned by our company and its affiliates. In case of any inconsistency between the content of the instruction manual and the actual product, the actual product shall prevail. Our company reserves the right to interpret and modify the instruction manual within the scope permitted by law.

If you have any objections to any content or clause of the instruction manual, please submit a written objection to our company within seven days after purchasing the product. Otherwise, it shall be deemed that you have agreed to, understood and accepted all the contents of the instruction manual.

Please note that once the product is connected to the network and activated or the value of the product is significantly depreciated after trial. Consequently, once you have actually activated or tried the product, you will no longer be entitled to the seven-day unconditional return service. Connecting the product to the network or binding the Robint APP account are both considered product activation.

## Product List

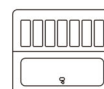

Main box x 1



Electrical Charge x 1



User Manual x 1

## Product specifications

Power supply: 12.00V—3.0A


Indicator light explanations:

- Yellow: Powered on but not connected to the network.
- Green: Powered on and connected to the network.
- Blue: Ejecting the medicine compartment.

## Operation Instructions

### ① Connect to Power

Connect the main box to the power. The indicator light turns on, system starts automatically.

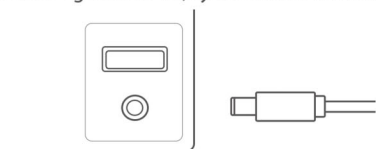

### ② Download Robint APP

Please scan the QR code on the below to download the Robint App, also can search "Robint" in App store.




### ③ Register your phone number and create a family

Sign up with your phone number, after that creating "Family".



### ④ Connect Internet

Please turn on Bluetooth, Search for the corresponding machine number, then select internet.

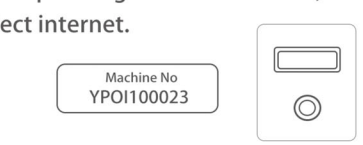

Note: If network configuration fails, repeat the above steps or contact customer service for consultation.

### ⑤ It can be configured in the app

Medication reminders can be created in the APP, and the medicine box can be controlled to eject the medicine compartment for placing medicines.



When power is off and the medicine compartment is removed. When power is off and the medicine compartment is removed.

## Warranty regulations

For non-human damage product performance failures within the warranty period, enjoy the following three guarantees services:

≤ 7days

Within 7 days of purchase, if there is a quality problem, consumers can choose to return or replace with the same model.

≤ 15days

Within 15 days of purchase, if there is a quality problem, consumers can choose to replace with the same model.

Note: During the warranty period, returned products will be inspected by authorized service providers. Defective parts will not be returned. Products without faults will be sent back as is. For online purchases, the purchase date is the product receipt date.

This certificate cannot be altered; otherwise, it is invalid.

## Non-warranty regulations

Situations not covered by the "three guarantees" (including but not limited to):

No "three guarantees" certificate or missing valid invoice or beyond the warranty period.

Damage caused by force majeure such as fire, flood, and lightning.

Damage caused by improper use, maintenance, or storage as per the user manual.

Damage caused by disassembly and repair by non-authorized service providers. Mismatch between product model/code on warranty certificate and the actual product.

Normal discoloration, wear, and consumption during use.

Inability to use due to consumer's own network issues. If the returned power supply has obvious hard object damage, cracks, broken pins, severe deformation, power cord damage, breaks, exposed cores, etc., the power supply is not covered by warranty or replacement.

In such cases, seek solutions from responsible parties, choose paid services, or purchase separately.

## Product three guarantees certificate:

Customer service name .....

Contact number .....

Address .....

Product model .....

Purchase date .....

#### FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

To maintain compliance with FCC's RF Exposure guidelines, This equipment should be installed and operated with minimum 20cm distance between the radiator and your body: Use only the supplied antenna.