


Wireless charging skin scrubber with mist

User Manual

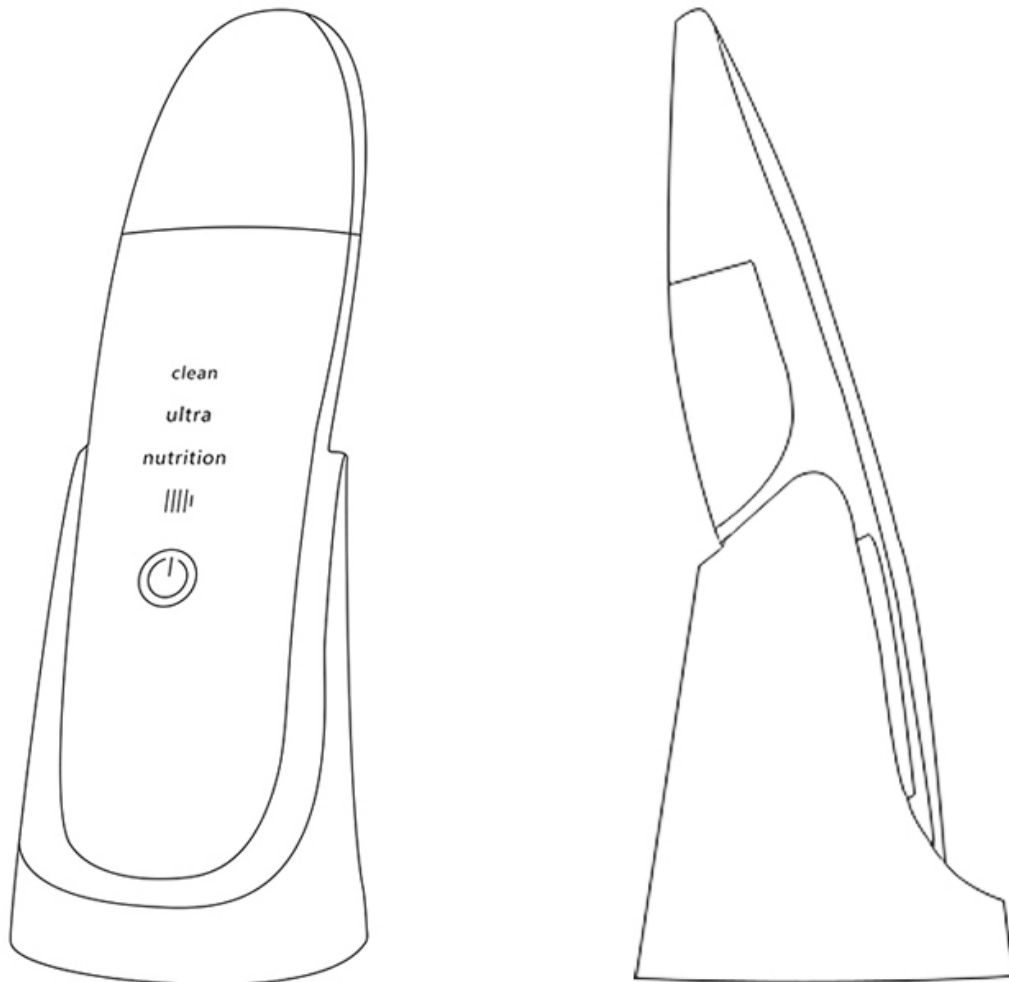

Attention

1. Do not exceed the scope of use specified in this user manual;
2. Pregnant women should follow doctor's advice
3. Do not use near area with metal prosthesis, such as eye balls,
inner ear, bone areas and joints;
4. Prohibit use in vein or arterial pathology areas or Patients with
Thrombophlebitis or Blood Coagulation Precedents;
5. Avoid physical impact;
6. Only use according to the user manual; Prohibit disassemble,
transform, maintenance or retrofit the product;
7. Please store the device in a clean and stable place: Do not store in
a place overheated or over cold; Keep away from fire, such as stove;
8. Charging base is not waterproof; Keep the charging base dry & clean;
9. Unplug the device after finish charging;
10. Keep the device away from children.

Introduction

Thank you very much for purchasing our products and your confidence in us; We hope our product can bring you a perfect experience.

Please read the user manual and all related information carefully and keep it for future reference.

Product Specification

N.W.:205G

Measurement: 21.5*9.8*8 cm

Working Environment: 0-35°C

Input: AC 100-240V, 50/60HZ; Output: DC-5V, 1A(Rechargeable)

Power: 2. 2W

Battery Volume: 800mAh

Ultrasonic Frequency: 24kHz

Blue Light Wavelength: 465±3nm

Waterproof Level: IPX5(Device only, charging base is not waterproof)

Accessories

1. Main Device
2. 5pin USB Charging Cable
3. Charging Base
4. User Manual
5. Gift Box
6. Water Transfer Bottle

Features

1. Standby Mode: Press the button once to turn on the device and enter STANDBY MODE, you can see the lights for "clean" , "ultra" & "nutrition" flashes.



2. Mist Mode: Press the button twice to enter this Mode, you can see the "clean" light flash and mist spray.



3. Continuous Wave Deep Cleansing Mode: Press the button third time to enter this mode, you can hear continuous sonic vibrations and see the mist spray, also the "clean" indication light always on.



4. High Level Deep Cleansing Mode: Press the button the forth time to enter this mode, you can hear continuous sonic vibrations and see "ultra" indication lights always on.



5. Sonophoresis Mode: Press the button the fifth time to enter this mode, you can hear continuous sonic vibration and feel ions. " nutrition" indication light always on.



6. Standby Mode again: Press the button the sixth time, enter Standby Mode again. No functions now. Press again functions circles.



7. Turn off the device: Long press the button for 1s to power off the device.



8. The device is set to work 10mins and then turn off.

How to charge

1. When low battery, the RED indication light lights up and device turn off automatically;



2. Put the device on the charging base and make sure well connected. Then charge it with USB cable. During charging, Red light flashes; When fully charged, light turns into Green and always on.



During charging, Red light flashes;



Finish charging, Green light always on.




Function Induction

1. Ultrasonic deep cleansing technology-The Ultrasonic Skin Scrubber
produces 24,000Hz per second of ultrasonic waves that effectively loosen clogged pores and exfoliates the skin to remove dead skin cells, painlessly.
2. Firms and lifts skin-micro-current pulse technology intermittent
ultrasonic helps boost cell collagen regeneration and activates cell vitality, ensuring your skin remains moisturized and healthy.
3. Advanced ion technology-The positive ion helps remove dirt and deeply cleanses the skin. It accelerates the speed of skin cell metabolism, softens old keratin, restores skin elasticity and helps keep skin feeling tight and smooth. The negative ion nutrition leading-in function helps break down and refine larger molecules such as collagen to massage and support your skin to easily absorb nutrients.
4. Bio Blue Light: wavelength $465\pm 3\text{nm}$: diminish inflammation, kill bacterium and propionibacterium acnes, remove acne and scars etc.

TIPS FOR EFFECTIVE CLEANSING

Deep Cleansing:



1. Apply skin cleansing serum/oil on your face, keep your skin always moist.
2. Turn on the device. Choose the first "Continuous Wave Deep Cleansing Mode" or second "High level Deep cleansing Mode", hold the unit with the control panel is away from the face. Cleanse the face according to the diagram on the following page. Treatment time is typically 3-5mins each time, no more 5 times per week. We recommend start with the first "Continuous Wave Deep Cleansing Mode" .
3. Wash your face with water after treatment.

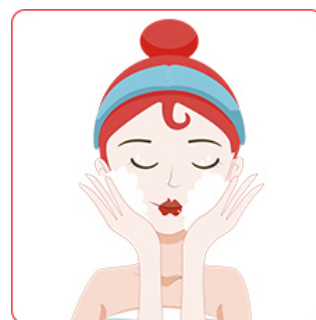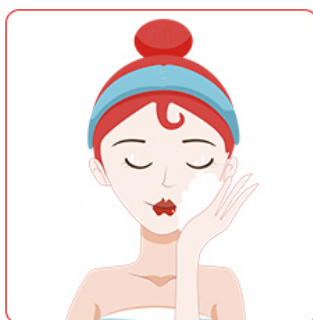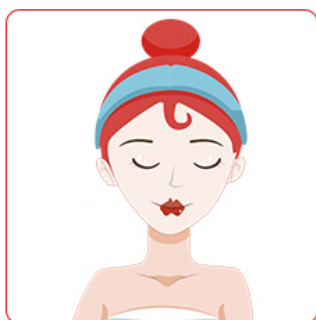

Cleaning step:

We recommend using the device according to this diagram, as it encourages the body's natural lymphatic drainage.

1. Sweep the machine from the neck upwards onto the face in horizontal strokes. Repeat on both sides of the face.
2. Methodically work over all areas of the face as shown in the illustration.
3. Change the direction of the handset and glide in vertical strokes starting from the neck area. Work from the jaw to the hairline.
4. Cleanse around the curve of the nostril using the curved inner edge of the blade.

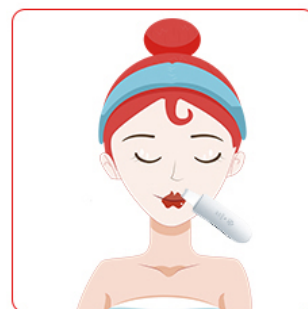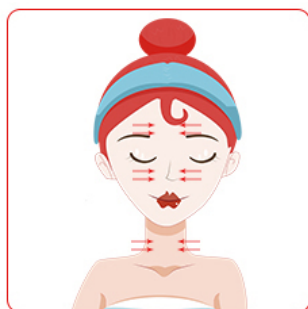

Use tips:


Use the scooped out middle section of the blade to cleanse the chin area under the lip with vertical strokes.

Select the first "Ultrasonic Deep Cleansing Mode" to work around the eyebrow area. Avoid the delicate under eye area when using the second "High Level Deep Cleansing Mode" .




Lift up & Sonophoresis

1. Apply your favorite serum or moisturizer on your face. Hold the device so that the blade is directly against the skin.
2. Turn on the device and the fifth "Sonophoresis Mode" .
3. Using gentle pressure, glide the blade over the surface of the skin using backwards movements. For best results, work slowly moving the blade at all times.
4. 3-5 minutes is enough for each treatment. Use for 3-5 times per week.

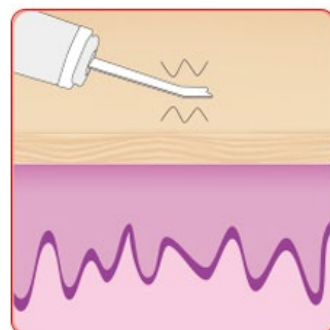


Tips For Lift up & Sonophoresis

Beginning on the neck area, with gentle pressure, glide the blade over the surface of the skin, using a backwards movement, working from the center of the face outwards.

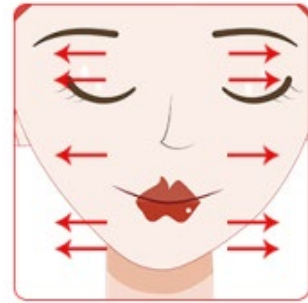

Continue upwards from the jawline to the hairline. Sweep over the face until all product is absorbed.

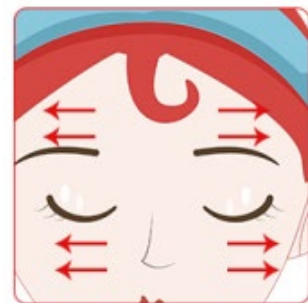

Maintenance


- 1 After each use, cleanse the device with a wet cloth and then keep it in a dry and cool place;
- 2 Before cleansing the device, make sure it is turned off and do not connect with power supply;
- 3 Personal use only; Do not share the device with others;
- 4 Keep this User Manual for future reference.

Failure

If the device do not work normally, please check following:


- 1 Make sure the device fully charged;
- 2 Make sure press "Power" button and the working indication lights lights up;
- 3 If no reaction after turn on the device, maybe low battery, please recharge the device and then try again;
- 4 If no action for 30s, the device will turn off Automatically to protect circuit.

How to add water

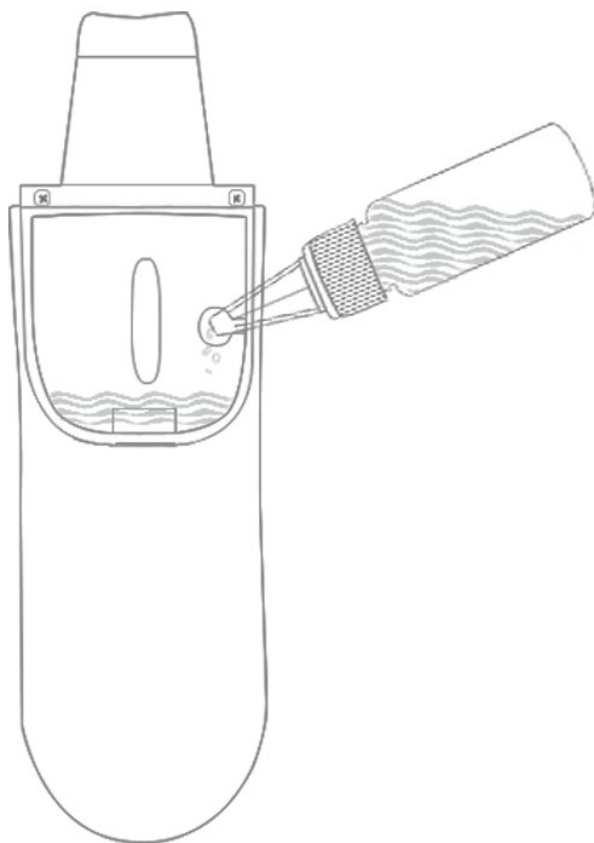

Waterproof

Waterproof: IPX5

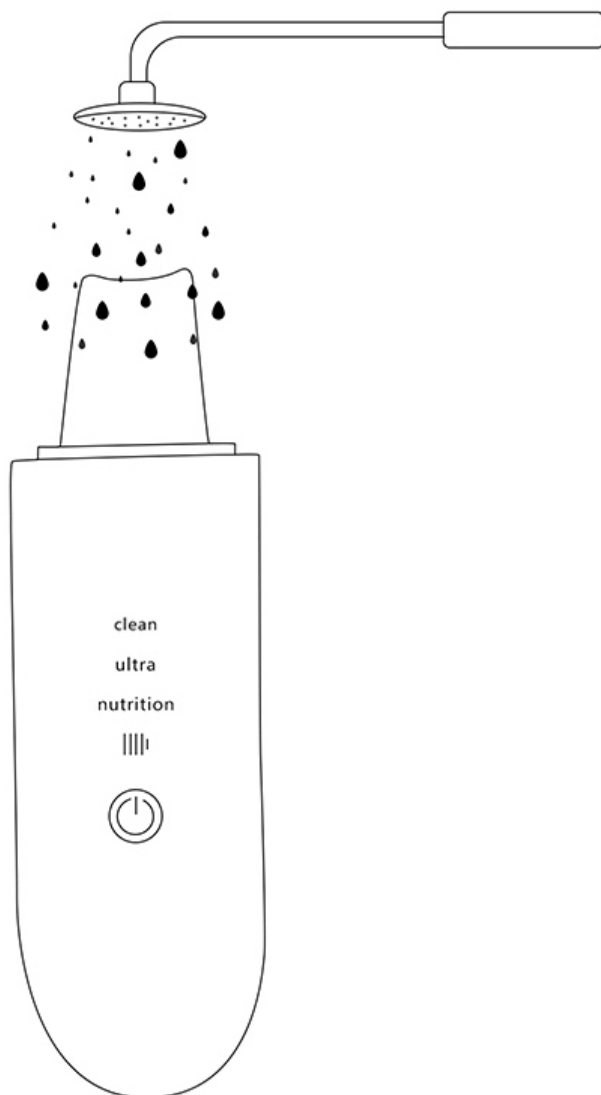



This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and

(2) This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- - Reorient or relocate the receiving antenna.
- - Increase the separation between the equipment and receiver.
- - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- - Consult the dealer or an experienced radio/TV technician for help.

This equipment complies with the FCC RF radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20cm between the radiator and any part of your body.