

Precautions

1.If the bluetooth is unable to connect to the Bluetooth® MiNi Massager for a long time, please close the Bluetooth® MiNi Massager software, then restart the software after 5s. Thesystem will search the Bluetooth® MiNi Massager and connect automatically.

2. When there is a phone call or message to come in, the phone automatically exit with operating interface, while the host is stillrunning. After answer the phone, it' s will return the operating interface, and massage to continue operating.

3.During the on-line operation, if the software has failed to search the host within 20s, the software will prompt that "It can't find the device. Please make sure the host power supply has been opened, and search again equipment".

4. Waiting for 30s to restart the host machine after power off,which makes mode reested easily.

5.When turned on the host, if the host don' t received the phone Bluetooth signal for a short time, or set aside after use, the host will power off automatically after five minutes later.

6.Host indicators:

	External USB Power				
Red	Light Battery charging	Off Battery charged	off	off	
Blue	off	off	Light Standby	Off Off state	Flicker Massaging

Figure IV

Use instructions

Introducing the distribution of the human body acupuncture point to the user in the guidance of the software, and the user can refer to the corresponding acupoints pasteur patch, in order to achieve the best effect.

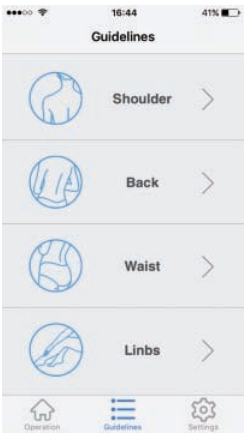


Figure V

Setting interface description

In the setup interface, the user can set the current name of the host device so that the user to distinguish between different Bluetooth® MiNi Massager. However, you need to re-pair connection to the host device after changing the name.

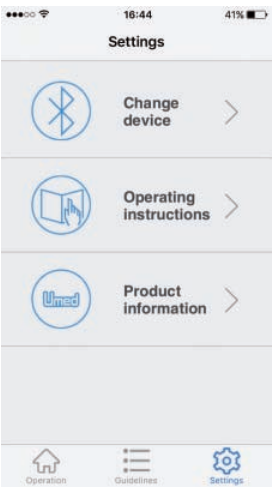


Figure VI

Charging method

Using the USB cable to connect the host and the computer, it is charging when the host red light, and the charging is complete when the red light off.

Tips

1. Prohibit the use of this product when using a pacemaker or other medical electronic equipment to continue the life;

2.Disable risk groups;

3. Can not be used on a wound, the body weak parts or active skin infections;

4. Female menstruation, pregnant women, hypersensitivity to the product, should be caution to use the product;

5. Prohibit two mutually adhesive patch when using the Bluetooth® MiNi Massager;

6.Relax while charging is not recommended;

7.Do not use this product when you are on the plane;

8.Do not connect USB cable charging when the computer starts;

Maintenance

1. When the host isn't used for a long time, it is recommended that charger the battery and check the host regularly, in order to prolong the service life of the battery;

2.Pay attention to Anti fall, waterproof, Anti high temperature, anti-corrosion;

3.Keep the product clean, it' s can use more soft cloth to wipe the host;

4.The dirt on the patch are available with the clean water.

Fault and Maintenance

Major failures	Possibilities	Solution	Remark
The light is off after the host machine is turned on	Battery exhausted	Charge the battery	
The connection is failed between the phone and the host machine	The equipment is out of the operating distance (> 6 m) Bluetooth signals are masked	According to Bluetooth, it can short the operating distance It's should avoid large obstacles to blocking	
	Low battery of the host machine	Charge the battery of the host machine	
The connection between the phone and the host machine is successful, but the host machine does not work.	The connection of the patch and the host machine is incorrect There isn't viscosity and conductivity for the paster used for a long time.	Reposted the patch in a right way Replace the paster	If the problem can not be solved, please call thecustomer service hotline of the Umed: 0755-26729100

Cautioned

1. Please strictly follow the Operating Instructions to avoid the adverse consequencesof improper use;

2.Due to personal physical differences, please adjust the position according to the feeling of the massage, and avoid to cause discomfort by high intensity.

Please keep this manual for inspection at any time, and read carefully before using the product.


Warranty card

Warranty Provisions
1.The warranty period is one year, effectively from the date of purchase.
2. The warranty scope only refers to the problems of Umed branded products caused due to the production, techniqueand quality.
3. The cases of the following circumstances, even during the warranty period, should also be repaired by pay.
a. The damage caused by usage, maintaining and storage notaccording to user manual.
b. The damage caused due to incorrect operation, such as power is not normal or broken(eg: drop)
c. The damage caused by dismantling of maintenance people who are not appointed by Umed brand.
d. Those who cannot present this warranty card and a valid purchase invoice.
e. Those Warranty recorded Information not conform to the actual commodity remark, or altered identification happens.
f. Damage due to force majeure (such as lightning, earthquakes, fire, flood, etc.)
4.During warranty period, the manufacturer will replace the faulty part free, but the freight out and home shall be borne by the customer.
5. When contacting brand distributor or manufactures for warranty repair, please present this warranty card and purchase invoice.

Item No.	H-LF- II	Purchase Date	
Warranty Period	one year, effectively from the date of purchase		
	Fill In By Customer		
Information Of Customer	Name		
	Phone		
	Address		
	Fill In By Seller		
Fill In By Seller	How to buy	Entity Shop	Online Shop
	Shop name		
	Contact Person		


Qualified Certificate

Inspector: _____



This product is strictly guaranteed according to national Three Guarantee regulations.

MINI MASSAGER



4.0
H-LF-II

MINI massager

USER MANUAL

Operating Principle

Umed Bluetooth® MiNi Massager adopting international leading technology, Send commands through the mobile phone Bluetooth module to control the host for a deep massage. The device adopt international advanced micro control chip technology and bluetooth technology, produce electronic pulse wave to Stimulate muscle, which make people to feel the feeling of tappings kneading, poundings tuggings pou&mox complex , Effectively relieve fatigue, alleviate muscle pain and other symptoms. Anytime, anywhere, let you relax.

Product Specifications


Host machine power	3.7V/180mA rechargeable Lithium battery
	Sleep Current < 100uA
	Working current < 60mA
	Standby Current < 12uA
Charging current < 120uA	
Product name	Bluetooth® MiNi Massager
Product model	H-LF-II
Phone requirement	Android 4.3 or more systems, IOS7.1 above, Bluetooth function
Operational mode	6 modes (tap+ knead+ pound+ tug+ pou&mox+ complex)
Level	10 levels adjustable
Output frequency	< 100Hz
Charging time	Charging 2.5 hours, continue working more than 20 hours
Time stalls	3 modes (15 minutes,30minutes, no time limit)
Size	58*58*14mm
Weight	20g
Implementation of standards	GB4706.1-2005/GB 4706.10-2008

Configuration list


Name	Number
Host machine	Standard
Paster	Standard
USB charging cable	Standard
Bluetooth® MiNi Massager software	Standard
Manual	1
Warranty Card	1
Certificate of compliance	1

Product configuration and function


The Bluetooth® MiNi Massager is consists of host machine, paster,USB charging cable and Bluetooth devices (Android phone orIphone). As follow:



1. Host machine
The host machine and paster are the main components of the Bluetooth® MiNi Massager.Power switch of the host machine is used to control the machine power.



2. Paster
Natural materials, which is installed on the host and affixed to place where you are uncomfortable.



3. USB charging cable
Connect to the host machine and the computer for charging.

FCC CAUTION

§ 15.19 Labelling requirements.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.§ 15.21 Information to user.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

§ 15.105 Information to the user.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

-Reorient or relocate the receiving antenna.


-Increase the separation between the equipment and receiver.

-Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

-Consult the dealer or an experienced radio/TV technician for help.

Instructions

1.Software installation instructions Scan QR code and install software



Scan QR code and download the software

2. User interface description on the phone

A. Open the box, install paster to the device, and then put the paster on the body which is unwell.

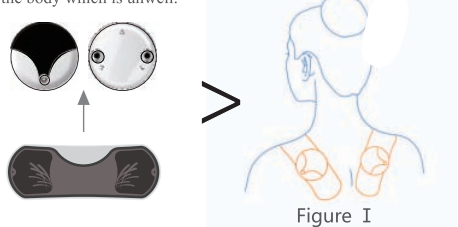


Figure I

B. Press the host machine power switch. (The button of the host machine is use the touch switch, press for power, successive double-click for power-off).

C. Open the bluetooth function in the "Settings" , click the icon on the desktop. Entering the software interface, the system search the bluetooth portable devices automatically. As shown in figure 2.

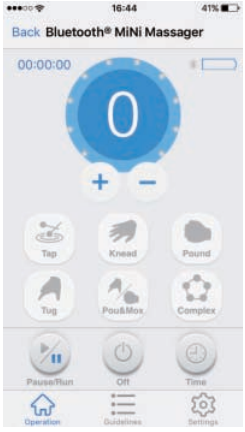


Figure II

D.Enter the main operation interface, select the corresponding mode, Press the button of increase the output intensity, begin to experience.

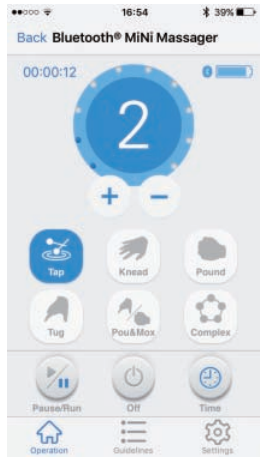


Figure III

There are "off" button、Start/pause button、MiniShutoff button、mode button, "+" button, "-" button on the operation interface;

"off" button——shutdown device , When using again, need to press the switch on the device manually;

Pause button——Click once, the device suspend the output, if you click on any pattern or gear switch, you can' t get any reaction, because the function is not effective; When clicked again, the device regain model and gear continue to work;

Model button——Click on the mode button, device can switch to the corresponding mode: tap、knead、pound、tug、pou&mox、complex;

Strengthen button -- click once, output can strengthen 1, until to 10th;

Weaken button——click once, output can weaken 1, until to 0;

E.Choose the appropriate mode, press function adjustment button, the weakest strength is 0, users can adjust to the appropriate state;

F.After finish relaxing, please click on the "off" button on the interface, or closed the device switch, and then remove the paster and separated from the device, to prevent dust damp cover paster with protective film.

G.Closed the device switch, and then remove the paster and separated from the device, to prevent dust damp cover paster with protective film.

Umed®

Company:Shenzhen Youyeah Health Technology Co.,Ltd
Add: NO1.1/F.2th building,Pingshan Minqi Science park,
Taoyuan Rd.Nanshan,District,Shenzhen,P.R.China
Customer Service Hotline:0755-26729100
Umed's website:www.umed.com

1

2

3

4