

ZAYATA

Smart Pill Dispenser


Quick Start Guide

Model:M126




■ Parts & Controls

- 1.Top Cover
- 2.Cover Handle
- 3.Dispense Door
- 4.Dosage Window
- 5.Light Indicator
- 6.Magnet
- 7.Dosage Ring
- 8.LCD Display
- 9.Button A
- 10.Button B
- 11.Button C
- 12.Removable Tray
- 13.Cover Lock
- 14.Battery Door
- 15.Battery Door Latch
- 16.Speaker

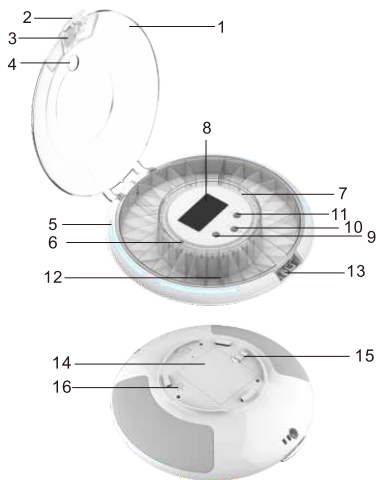

■ LCD Display Instruction

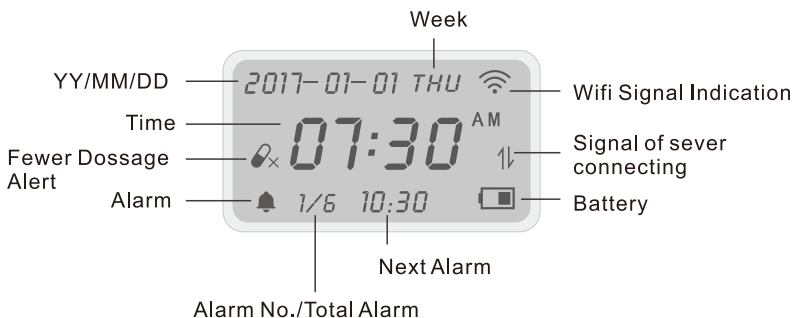

■ Inserting Battery

1. Slide Battery Door Latch, then remove Battery Door.
2. Insert battery with correct positive and negative poles




■ Open or Lock Top Cover

1. If locked, locate Cover Lock on bottom of unit and insert the included key, then turn clockwise to unlock.
2. Put Pill Dispenser on table, press Cover Handle to release the latch, then lift to open Top Cover.
3. To lock Top Cover, please close Top Cover first, press Cover Handle to get latched, then turn counterclockwise to lock with inserting the key from bottom side.



■ Place the Dosage Ring

1. Provided dosage rings with the Pill Dispenser for 1-6 doses per day.
2. Set correct Dosage Ring on center of tray, and fold the tab down to secure in place.
3. Place pills in correct dosages in each compartment.
* Do not place pills in compartments corresponding to red markings on Dosage Ring.

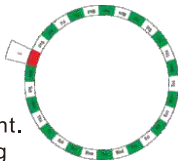

■ APP Download & Smartphone Connecting


1. Download from link <http://qr.zayata.com/qr/>. Or scan below QR code.
2. Open the APP, log in with ID No. & Password. If you are a new user, please sign up first, then log in.
3. A scan interface show up at first time after log in, please scan the QR code at the back of device to set up bind.





Download the QR code






Scan the QR code with smartphone

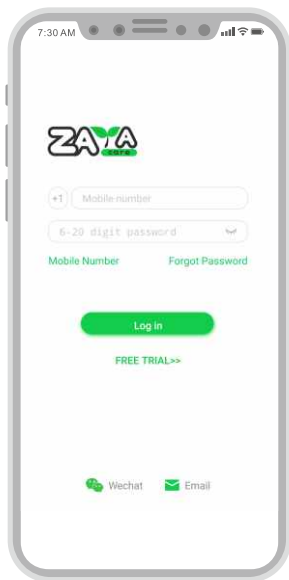

4. WiFi Connecting

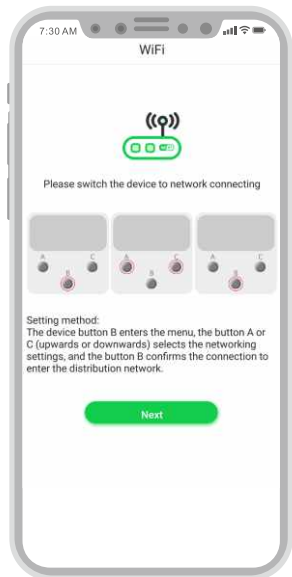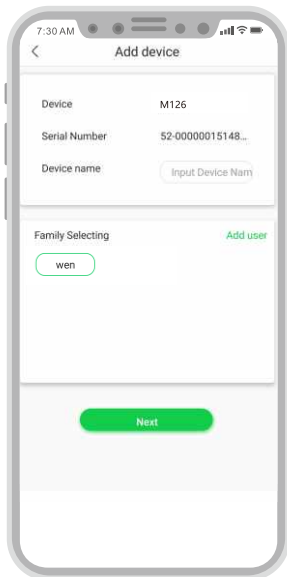


5. Please follow the stand-alone device steps to set the device connecting with Wifi.

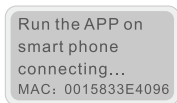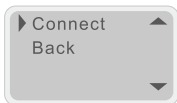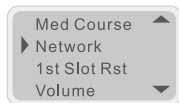

6. Set up WiFi of device to WiFi Router

1. Make sure smart phone already connected with WIFI router (which can connect to Internet).
2. The APP will get SSID of WiFi automatically from smart phone.
3. Input same "WiFi password" as smart phone.

- 4.Waiting for smart phone and device synchronization result.
If connecting overtime, try again to reconnect.

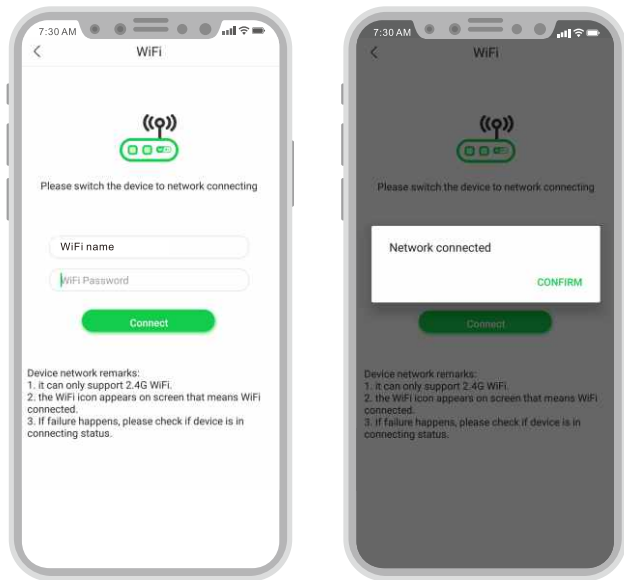

7. “Successful” will show up in device and APP at the same time.
It means the device will connect WiFi automatically.
- 8.If the connecting is overtime in device or APP, it will show you
“Connection Failed”. Please try again to reconnect.

Successful!


Faile!

9. Open the Top Cover, take out the Removable Tray(11). Set the red mark of Dosage Ring in first grid (no compartment), then fold the tab down to secure in place.



10. Set up new Wi-Fi connection for device. Please follow below process for setting, Click "Device", choose "Setting", and then "WiFi".

11. The sequence of placing pills inside.

Please pay attention that place pills with counterclockwise direction.

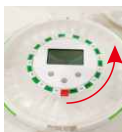



12. Open the APP to set up medication plan after placing the pills finished.

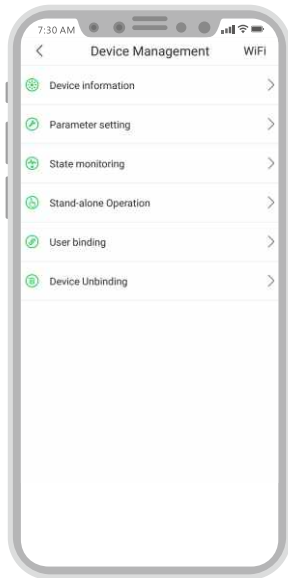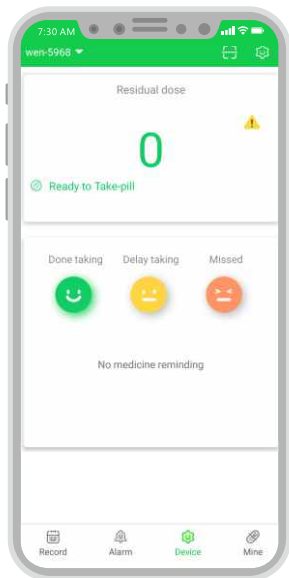

13. Set up Medication Plan


In Alarm interface, click “Edit” to enter alarm settings, click “open alarm” to open, set up the time as needed and click “confirm” to save the alarms.

“Alarm Expire Time” means effective period of alarm time. If you set up from 2017/09/27 to 2017/11/27, the alarm will ring during this period. It will not ring out of this period.

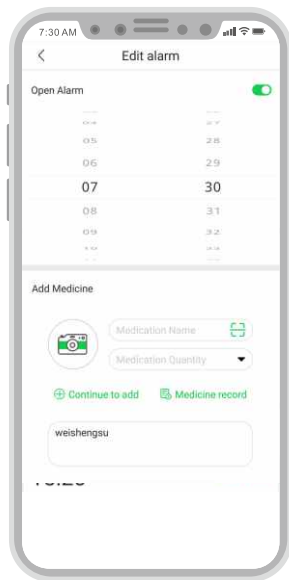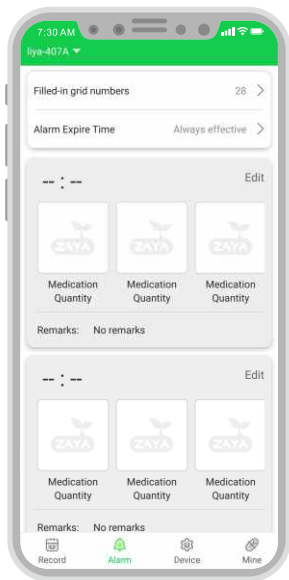


14.1 to 6 medication alarms, setting time according to your own needs.

If need to close the alarms, click “Open Alarm” to close. Then click “Confirm” to save the settings.

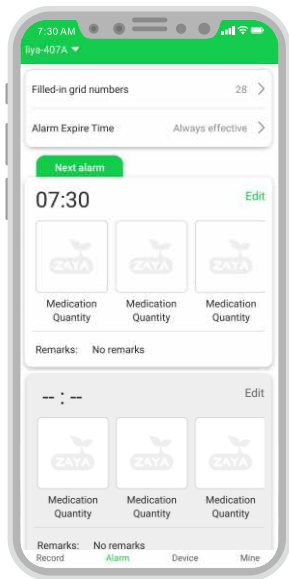


14.Remove binding

Click “Device”,and then “setting”,then “Remove binding” .
When remove binding,the parameter will reset at the same time.

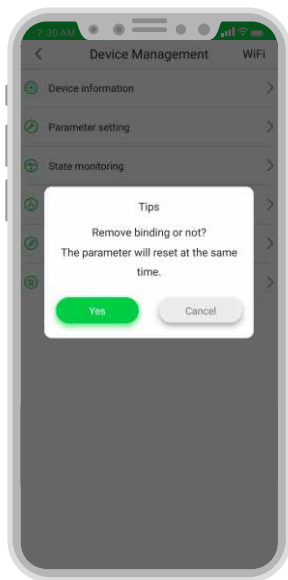


■ Stand-alone Operation Instruction

1.Alarm Setting

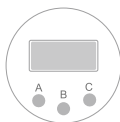

Button A means Up. Button B means Confirmation. Button C means Down.

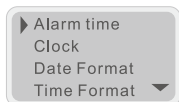


Press Button B to enter Menus. Choose “Alarm Time” for alarm setting. Then press Button B, the hour will flash, press Button A or C to adjust. After you setting up the hour, press B to go on minutes, press Button A or C to adjust. Press Button B again means 1st alarm setting “Done”.



Setting up other alarms, please follow same step above.

Press Button A or C to “Return” tip, then press Button B to get back to Menus.


2.Clock Setting

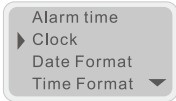


Press Button B to enter Menus. Choose “Clock” for setting.

Press Button B When the digit is flashing, press Button A/C to adjust.

Press Button B to adjust next digits, save and exit.

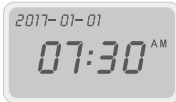

3.Date Format Setting

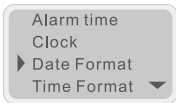


Press Button B to enter Menus. Choose “Date Format” for setting. When the digit is flashing, press Button A/C to adjust. Press Button B to adjust next digits, save and exit.

4.Time Format Setting



Press Button B to enter Menus. Choose “Time Format” for setting. When the digit is flashing, press Button A/C to adjust. Press Button B to adjust next digits, save and exit.

5.Med Course Setting

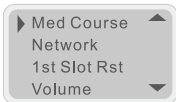


Press Button B to enter Menus. Choose “Med Course” for setting. When the digit is flashing, press Button A/C to adjust. Press Button B to adjust next digits, save and exit. The function is for limited alarm date setting from start to end.

6.Network Settg




Press Button B to enter Menus. Choose “Network” for setting. Press Button A/C to choose “Connect”, press Button B for WiFi connecting.

7.1st Slot Restoration

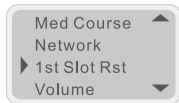


Press Button B to enter Menus. Press Button A/C to choose “Confirm”, then press Button B. Then the device will automatic restoration. Or Long press Button A/C for more than 5 seconds to 1st Slot RST.

8. Volume Setting

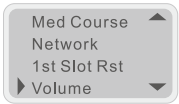


Press Button B to enter Menus. Press Button A/C to Volume, Press Button B to choose. Press Button A/C to choose the volume you want. Then press Button B to save and exit.

9. Language Setting

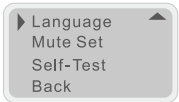

Press Button B to enter Menus. Press Button A/C to “Language”, then Press Button B to choose and exit.

10. Mute test

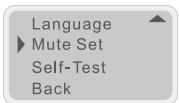

Press Button A/C to “Mute Set”, press Button B to choose. Then press Button A/C to choose the setting you want, press Button B to confirm.


1. Non-interruption ON-OFF. Press Button B to choose on or off non-interruption mode.

2. Start time Setting. Move cursor to "Start time" and press Button B to set the time, same as alarm setting.

3. End time setting. Move cursor to "End time" and press Button B to set the time, same as alarm setting.


4. Exiting setting. Press Button A/C to "Back", Press Button B to choose.

11. Self-Test

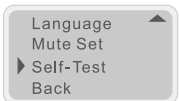

Self-test is only for factory testing.

12. Back

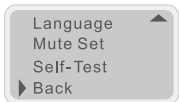

Press Button B to enter menus. Press Button A/C to adjust to “Back”. Press Button B to get back.

■ FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Warning:

§ 15.21 Information to user.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

§ 15.105 Information to the user.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

* RF warning for Portable device:

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment.

This equipment should be installed and operated with minimum distance 20cm between the radiator & your body.