



Smart Moxibustion Device

Product Manual

1. Product list

- Moxibustion device
- Moxa filler
- Charging base
- Charging wire
- Product Manual
- Certificate
- Warranty card

2. Operation steps

- 1) Remove the device from its charging base

Button



Function of the button


- Turn on
- Temperature adjustment
- Turn off

2 Load the moxa filler

Take out the moxa filler



Tear off the release paper on the ring surface of the moxa filler, clinging to the bottom of the device.

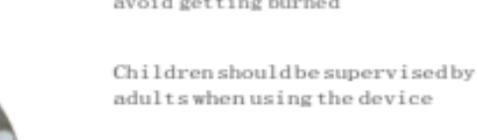

Ring surface of the moxa filler



3 To turn on the device

Press and hold the button for 2 seconds

Take out the moxa filler



Patch the moxa filler on the area you want the device to work on to avoid getting burned



Children should be supervised by adults when using the device

Children should not be allowed to play with device to avoid getting burnt.



4 Find the acupoint

Press and hold the button for 2 seconds until the light indicator is off.

Charging probe of the device.



Conductive needle of the charging base.



Children should be supervised by adults when using the device

Children should not be allowed to play with device to avoid getting burnt.



5 There are two ways which you can use to adjust the temperature.

Press the button once for gentle heat

Pressing the button twice for warm heat

Pressing the button for the third time for hot heat

Options 2

Download iJoou from Appstore, or Android

<https://intl.ijoo.com/download.html>



1) Activate Bluetooth, please charge Moxibustion device firstly to activate Bluetooth, as Bluetooth was set to hibernate originally

2) Turn on Bluetooth and connect device to app

3) Use the app to select the perfect temperature

6 To turn off the device

Press and hold the button for 2 seconds until the light indicator is off.

Charging probe of the device.




Conductive needle of the charging base.



Children should be supervised by adults when using the device

Children should not be allowed to play with device to avoid getting burnt.

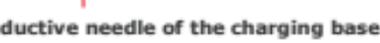

3. Charging indication

2) Plug-in the power cord and turn on the socket.

3) When the temperature light turns red the device has started charging.

4) When the temperature light turns green, the device is fully charged and will disappear automatically after one minute and it's ready to use.

5) If the temperature light is off or the red light twinkles after pressing the button, it indicates that the power capacity is low.

4. Important safety cautions!

1) Consult with your physician before beginning this program especially if you are suffering from:

hypertension, diabetes, cardiovascular, cerebrovascular disease.

2) When conducting moxibustion start with low temperatures and increase gradually to avoid discomfort caused by too much heat.

3) You should not conduct Moxibustion frequently at least once or two times a day.

4) Use the product with caution, if you feel dizzy, short of breath, starvation, thirst, fright or anger.

5) Do not conduct Moxibustion for a long period of time, take a break at an interval of 30 minutes.

6) Persistent Moxa therapy would have a better effect.

After Moxibustion, peel off the flower surface of the moxa off your skin and tear the ring surface of the moxibustion device slightly.

Smart Moxibustion Device

Warranty Card

M10, Central Park of High-tech Park in Nanshan District (No. 4 Kewen Road), Shenzhen

400-0755-929

Smart Moxibustion Device

Product Certification

Product name Smart Moxibustion Device

Product model ijoo-B02

Manufacturing Date _____

Inspector _____

Shenzhen Baichuan Technology Co., Ltd

FCC Warning Statement

Changes or modifications not expressly approved by the party responsible for compliance could void the user's ability to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

--Reorient or relocate the receiving antenna.

--Increase the separation between the equipment and receiver.

--Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

--Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and

(2) This device must accept any interference received, including interference that may cause undesired operation.

Maintenance date

Maintenance record

The following conditions are excluded from the warranty scope:

1. Products purchased from non-authorized ijoo online store or dealer;

2. Damage or fault caused by misuse, attempts to repair, or dropping;

3. Damage or fault caused by inappropriate operation, going against the manual;

4. Damage or fault caused by inappropriate/improper maintenance or storage;

5. Products repaired by unauthorized ijoo service center

Smart Moxibustion Device

Product Certification

Product name Smart Moxibustion Device

Product model ijoo-B02

Manufacturing Date _____

Inspector _____

Shenzhen Baichuan Technology Co., Ltd

FCC Warning Statement

Changes or modifications not expressly approved by the party responsible for compliance could void the user's ability to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

--Reorient or relocate the receiving antenna.

--Increase the separation between the equipment and receiver.

--Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

--Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and

(2) This device must accept any interference received, including interference that may cause undesired operation.