


# ***TOLOCO***<sup>®</sup>

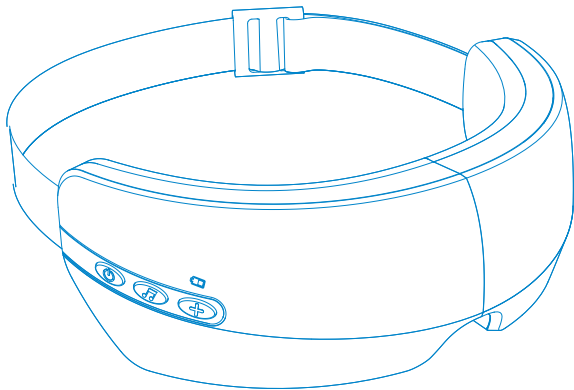

TL-EM08

**Eye Massager**

INSTRUCTIONS

### **Attention:**

**Please read and comprehend all the contents of this manual and operate the product according to the instructions.**


**DO NOT** use this product if have eye diseases like Glaucoma, Cataract, and Retinal detachment. If eye ever carried out surgery, it is prohibited to use this product.

- **DO NOT** clean the product with water, clean it with cotton soaked in alcohol.
- **DO NOT** disassemble the product.
- This product has a hot surface, people sensitive to heat should take extra care when using this.
- This product is not intended to be used by people who are physically weak, slow in response, or mentally disabled, unless it is used safely under the guidance or help of the person responsible for their safety. Unattended young children and disabled persons are prohibited from using this appliance, and children are prohibited from playing with this appliance.
- This product is intended for personal care use, it is not intended to replace medical care and medical devices.
- **DO NOT** disassemble the product. **DO NOT** impact or squeeze its interior lithium battery or expose it to fire. Do not use when severe swelling occurs.
- **DO NOT** use if the battery has been immersed in water.
- Lithium battery set inside the product. **DO NOT** place in high temperature environment.


- Consult your doctor if you have sensitive skin, allergic dermatitis, or orthopedic bio materials on face.
- Remove contact lenses, hairpin and hair accessories before using.
- Start the eye massager after putting it on, in order to avoid damage to the air bag due to over pressure.
- **DO NOT** use in high temperatures or humid environments. Avoid direct exposure to the sun and gas corrosion.
- This product will block the view. **DO NOT** use when walking or driving.
- Only use the charger (5V/1A) to recharge the eye massager.
- You must have your eyes closed while massaging.
- After massaging, the eyes are relatively blurred, and slowly become clear after closing the eyes for 1 minute. Run time: around 2-3 hours
- **DO NOT** use sharp objects to scratch this product.

## Product Structure

### Wireless BT Main Engine

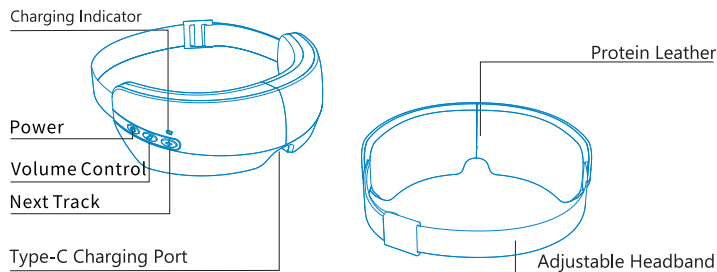


## Accessories

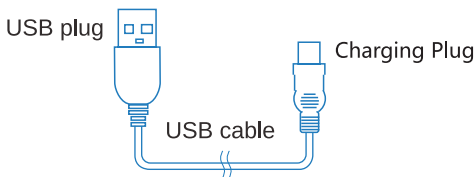

### *Preparation before use --- Charging*

- Before using, use 5V mobile charger for charging eye massager.
- Insert USB plug into the socket of the mobile charger.
- Charging Plug into the socket of the device, make sure the Charging Plug is inserted in the right direction.
- While charging, the red LED indicator will display. When done charging, the LED will turn green. The device will not work while charging is in progress.
- The device will not work while charging is in progress.
- The charge duration will be 2-3 hours. Run time: around 2-3 hours  
Disconnect the power adapter after charging is complete.

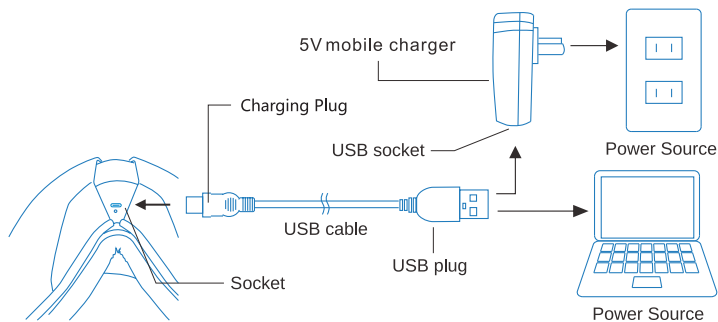


## Using Instructions

### Wearing the Device

Align the cloth cover to your eyes, put on the device, adjust the elastic band, adjust the device to the most comfortable position.

Hold the power button for 3-5 seconds to turn on the eye massager.

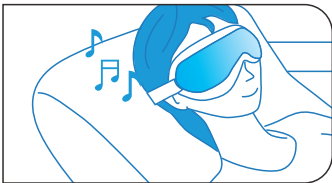

### ⚠ Attention:

Remove any glasses, contact lenses, hairpins and hair accessories before use to avoid discomfort. Do not operate the device with wet hands. Do not drop or step on the product.

### Key Control

[⏻] Long press 3 seconds to turns on/off Eye Massager.  
Short press switch massage mode

[🎵] Short press: next track

[+] Short Press: Volume adjustment (level 4 default-5-6 maximum-1-2-3). Long Press Turns Off music

### Wireless BT Music

BT pairing name: **Eye Massager**

- Wireless BT is **ON** by default.
- Massager will play built in music after start up.

- It will switch to Wireless BT device to play the BT music after mobile phone connection to the BT.
- It will automatically restore the built in music play after Wireless BT disconnection.

**Note:** Over time, repeated exposure to loud noise and music can cause hearing loss.


## Massager Programs

### Default Setting:

**Soft Mode:** Soft Air Pressure + Hot Compress + Music (15 minutes)

**Press POWER button to cyclely changing massager mode:**

Everytime short press/touch **POWER** button, Audio voice will confirm mode change. (music will stop during this operation, it will auto restore music after 2 senconds)

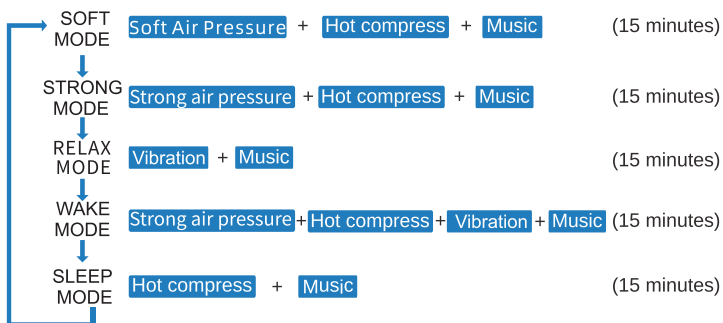

## Products Features



Hot Compress



Automatic Timing



Air Pressure



Vibration



Intelligent  
Voice Reminder



BT



180° foldable

## Transport and Storage

**Temperature requirement:** -40 ~ 50°C (-40°F ~ 122°F)

**Relative humidity:** ≥80%, No corrosive gas, good ventilation.

This product is equipped with simple anti-vibration facilities, suitable for air, railway, highway and ship transportation.

Rain, snow, upside-down and collision should be avoided.

When the storage period is more than three months, you should open the box for inspection, and do the power-on processing for half an hour.

## FAQ and Troubleshooting

Problem	Possible causes	Actions
There is a sound mechanical or air sound coming from device	Normal sound of air pump inflating, motor working, or air pressure	No need to do anything
Device turns off automatically. Device does not turn on. Buttons are not responsive. Air pressure is weak.	Battery is low	Charge the device

- If there are other problems, please cut off the power supply for the abnormality after being treated as above and then contact with the distributor or our company.
- Do not disassemble the product for maintenance by yourself.  
Our company shall not be liable for product damages or personal injuries therefore.

## Cleaning and Maintenance

### ⚠ Attention:

Turn off the power when cleaning, do not use corrosive detergents such as benzene or diluents to clean the product. If long time no use the device, please use your own anti-dust cloth to protect the device.

## Guarantee

This product comes with one-year warranty from the date of purchase.

### **Warranty will be voided due to any of the following:**

- Damaged due to the improper usage or storage.
- Damaged by attempting to repair or disassemble the device.
- Period of warranty exceeded.



## ***Products Features***

- Compact & Portable Facial Massage Therapy System
- Selectable Preset Eye Massage Modes
- Assists in the Relief of Eye Strain, Dry Eyes, Puffiness, Sinus Pressure, Headaches, Stress Relief and Helps to Sleep
- Acupressure: Simulated Pressure Pulsing
- Heating Mode: Hot Compress for Improved Blood Circulation, Helps Alleviate Dark Circles
- Pneumatic Mode: Simulates Pressure of Soft Hand Massages
- Music: Wireless Eye Massager with Soothing Sounds
- With Wireless BT Speaker to Play your Own Soothing Sounds
- Sleek, Next-Gen Fashionable Goggle Visor Style
- Simple One-Button Activation
- Built-in Rechargeable Battery
- Ultra-Quiet Operation
- With USB Cable
- Adjustable Headband Wrap Design
- Universal Size for All Ages: Adjustable Elastic Band
- Perfect for Use After TV & Computer Screen Viewing, Before Bedtime

### ***Wireless BT Connectivity:***

- Built-in BT for Wireless Streaming
- Works with Smartphone Only
- Wireless BT Network Name: 'Eye Massager'
- Wireless BT Wireless Range: 98'+ Ft.
- Wireless BT Version: 5.0

### ***What's in the Box:***

- Digital Eye Massager
- User Manual
- USB Charging Cable
- Desiccant

### ***Technical Specs:***

- Construction Material: Engineered ABS, Soft PU
- Recommended Massage Session Duration: 15 Minutes
- Eye Massage Pressure Strength: 85 kPa
- Rechargeable Battery: 3.7V Lithium, 1500 mAh
- Power Supply: 5V DC, 1A (USB Powered)
- USB Charging Cable Length: 1m
- Product Dimensions (L x W x H): 8.07in\*1.2in\*2.95in

## ***FCC Warning:***

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

### ***IC Warning:***

This device complies with Industry Canada's licence-exempt RSS standard(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference; and
- (2) This device must accept any interference, including interference that may cause undesired operation of the device. Cet appareil est conforme aux CNR exemptes de licence d'Industrie Canada.

Son fonctionnement est soumis aux deux conditions suivantes :

- ( 1 ) Ce dispositif ne peut causer d'interférences ; et
- ( 2 ) Ce dispositif doit accepter toute interférence , y compris les interférences qui peuvent causer un mauvais fonctionnement de l'appareil. This Class B digital apparatus complies with Canadian ICES-003. Cet appareil numérique de la classe B est conforme à la norme NMB-003 du Canada. The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.