


UREVO Care Wireless Relax Massager

User Guide

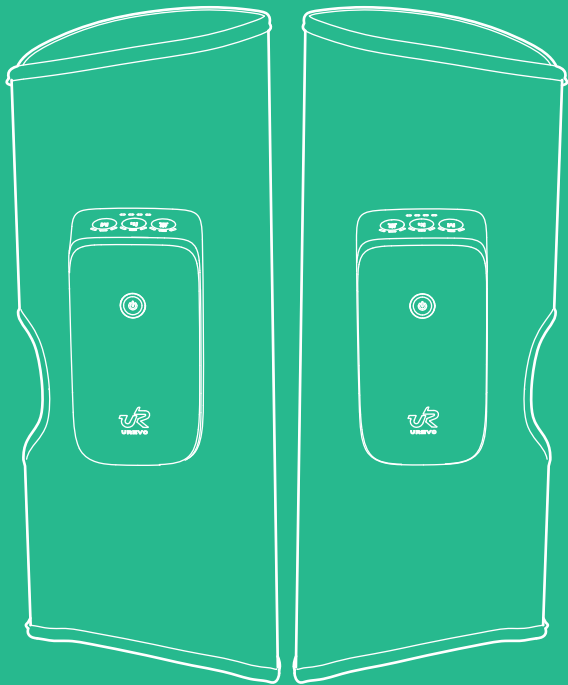

Scanning QR code
UREVO APP

US
Model. UCRM001



REGISTER NOW TO ACTIVATE 24-MONTH WARRANTY!

You'll receive a 24-month warranty and lifetime customer service after activation. This warranty applies only to eligible purchases from UREVO authorized distributors. For more details, visit urevo.com/warranty.

Customer Service

Email:

- **US:** services@urevo.com
- **EU:** service_eu@urevo.com (Response within 24 hours)

US Support Hotline: (1)844-988-2473

Hours of Operation: Monday - Friday (9 am - 12 pm; 1 pm - 6 pm PST)

IMPORTANT: Read all instructions carefully before using this product.

The specifications of this product may vary from this photo, subject to change without notice.

Please refer to the detailed manual on the **UREVO APP**.

Please open **UREVO APP** and go to "**Device details - Settings**" or "**Profile**" to view the detailed manual.

 **For more content, please visit : urevo.com**

IMPORTANT SAFETY INSTRUCTIONS

Read all instructions before using (this appliance).

DANGER – To reduce the risk of electric shock:

1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

1) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.

2) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

3) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.

4) Keep the cord away from heated surfaces.

5) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.

6) Never drop or insert any object into any opening.

7) Do not use outdoors.

8) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.

9) **WARNING:** Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

10) **WARNING:** Do not stand on or in appliance. Use only while seated.

11) This appliance contains batteries that are non-replaceable.

12) This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

SAVE THESE INSTRUCTIONS



IMPORTANT SAFETY INSTRUCTIONS

- 1.Do not use while it is charging
- 2.While using, turn off the power first if intending to stop using.
- 3.Do not use while driving
- 4.Do not use during sleep.
- 5.Only professional staff allowed to disassemble or repair this product.
- 6.Do not use when power cord is damaged
- 7.Do not use when water of other liquid leakage.
- 8.Do not use when the product is abnormal.
- 9.Stop using immediately if you have any discomfort.
- 10.Do not place this product in inflammable,explosive or damp environment.
- 11.Please avoid sharp instruments contacting with this product ,so as not to damage the product.
- 12.This product contains lithium battery,before abandoning the appliance,power off the appli-
ance and take out the lithium battery the battery must be disposed according to local regula-
tions,the battery can not be treated as domestic waste,do not throw the battery into fire.

 Warning of personal safety or property damage caused by incorrect operation.

NOTICE

IT SHOULD BE USED WITH CAUTION OR UNDER THE GUIDANCE OF MEDICAL DOCTORS IN CASE OF THE FOLLOWING SITUATIONS.

- Pregnant or postpartum women
- Suffering from bone disease
- Skin allergy (e.g.scald or burn)
- Patients with severe heart disease and abnormal blood pressure
- Patients with frost epilepsy
- The body has implanted devices(e.g.pacemakers)
- People with osteoporosis

Specifications

| | |
|------------------|------------------------------------|
| Product Name | UREVO Care Wireless Relax Massager |
| Product model | UCRM001 |
| Rated power | 5W |
| Rated voltage | 5V |
| Battery capacity | 2500mAh |
| Default timing | 20mins, heating 15mins |

FCC STATEMENT:

Operation is subject to the following two conditions:
(1) This device may not cause harmful interference, and
(2) This device must accept any interference received, including interference that may cause undesired operation.
(3) Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
Note:This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try the interference by one or more of the following measures:
to correct
• Reorient or relocate the receiving antenna.
• Increase the separation between the equipment and receiver.
• Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
• Consult the dealer or an experienced radio/TV technician for help.

This device complies with RSS247 of Industry Canada. Cetappareil se conforme à RSS247 de Canada d'Industrie. This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device. appareils radio exempts de licence. Son fonctionne-
mentest sujet aux deux conditions suivantes: (1) le dispositif ne doit pas produire de brouillage-
préjudiciable, et (2) cedispositif doit accepter tout brouillagerecu, y compris un brouillage
susceptible de provoquer un fonctionnementindésirable.
Radiation Exposure Statement: The product comply with the Canada portable RF exposure limit
set forth for an uncontrolled environment and are safe for intended operation as described in
this manual.

The further RF exposure reduction can be achieved if the product can be kept as far as possi-
ble from the user body or set the device to lower output power if such function is available.
Déclaration d'exposition aux radiations : Le produit est conforme aux limites d'exposition pour les
appareils portables RF pour les Etats-Unis et le Canada établies pour un environnement non
contrôlé.
Le produit est sûr pour un fonctionnement tel que décrit dans ce manuel. La réduction aux exposi-
tions RF peut être augmentée si l'appareil peut être conservé aussi loin que possible du corps de
l'utilisateur ou que le dispositif est réglé sur la puissance de sortie la plus faible si une telle fonction-
est disponible.
This Class B digital apparatus complies with Canadian ICES-003.
Cet appareil numérique de la classe B est conforme a la norme NMB-003 du Canada.