

# Minimalist MG04

MEDCURSOR

英文版说明书 成品尺寸: L85 x W85mm 展开尺寸: L85 x W170mm © copyright:zhao 20200828 改 材质: 白色道林纸



MEDCURSOR

## CONTENTS

Safety	03
--------	----

Product Care	10
--------------	----

Inside the Box	12
----------------	----

Product Features	13
------------------	----

Instructions Manual	16
---------------------	----


Troubleshooting	24
-----------------	----

Product Disposal Requirements	25
-------------------------------	----

Product Spec	29
--------------	----

EN


## Precautions

	This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance. User should safely operate the device and should know the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
	The appliance should only be used with the power supply unit provided with the appliance.
	This appliance contains batteries that are only replaceable by skilled persons.
	SRRC-certified model approval code, a.k.a. CMIIT ID, can be found on the product's nameplate.


3

EN


## Restrictions and Reminders

	Read the Instructions Manual carefully. Avoid long-time continuous use of this product to prevent muscle pain caused by excessive massage.
	Do not use this product while driving or conducting precision operations.
	Do not use on surgical wounds or damaged/broken skin or area with inflammation or congestion.
	Do not use if suffering from heart disease (or with a history of heart disease), acute disease, skin disease, malignant tumors, blood disease (such as bleeding), osteoporosis, soft tissue disease, stroke or other physical discomforts.
	Should not be used by anyone who have low tolerance to pressure, who feels discomfort or pain or anyone with physical discomfort or intellectual disability. Please consult your doctor if you have any medical condition that could be affected by the massager.


4

	Do not use if pregnant or when menstruating unless your doctor advised that it's safe.
	Not for children under 8. Children aged 8 to 12 are not recommended to use this product unless they have been given supervision or instruction concerning the use of this product by parents or guardians.

## Caution

	This is a household electronic healthcare product and shall not be used as medical device or for diagnosis, treatment or other commercial purposes.
	Before use, remove accessories around the massage area or other items that may affect the massage.
	Stop using if pain or discomfort occurs in the massage area.
	During use, if any issues occur, cut off power supply immediately and contact customer service. Do not disassemble or repair it by yourself.
	Please observe safety regulations.

## Product Care

	Clean the product with soft dry cloth. Do not rinse or wipe with water.
	Clean the surface of the product with soft wet cloth or a 75% alcohol pad. Do not rinse.
	The base of this product is not waterproof. Do not rinse.
	Charge the product with the charging device provided by our company. Practice safety precaution.
	Observe safety precaution to avoid electrical hazards or risks.
	Keep the product away from water, fire and corrosive environment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The statements should be displayed in the user manual:

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

## EN

### Must Know

This massage head is made of food-grade silicone/metal that can directly contact the skin. Please read the Instructions Manual carefully before use to avoid skin damage.

#### ※ Usage Frequency

Twice daily (morning and evening) for no longer than 10 minutes each is recommended.

#### ※ Cautions Before Use

1. Do not use with wet hands.
2. Stop use and seek medical attention if pain or burning occurs in the massage area.

#### ※ After Use

1. The massage area may become red and warm in a short time, which is normal and should recover soon.


### ※ How to Clean

1. Clean the product with soft dry cloth. Do not rinse or wipe with water.
2. Clean the surface of the product with soft wet cloth or a 75% alcohol pad. Do not rinse.

### ※ Others


This product may hum or make a low buzzing sound due to pressure or vibrating feature when it is being used. This is normal.

### PRODUCT CONFIGURATION

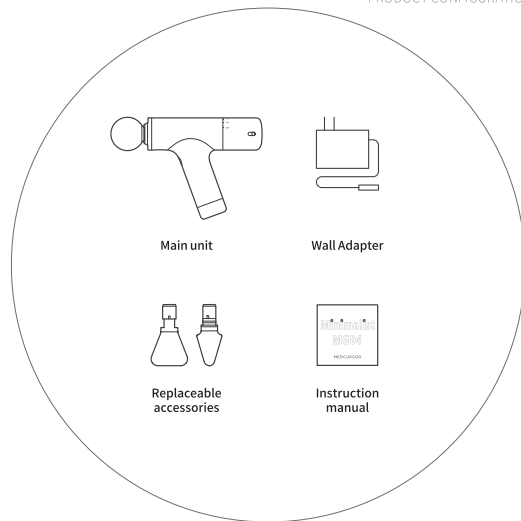




## PRODUCT FEATURES

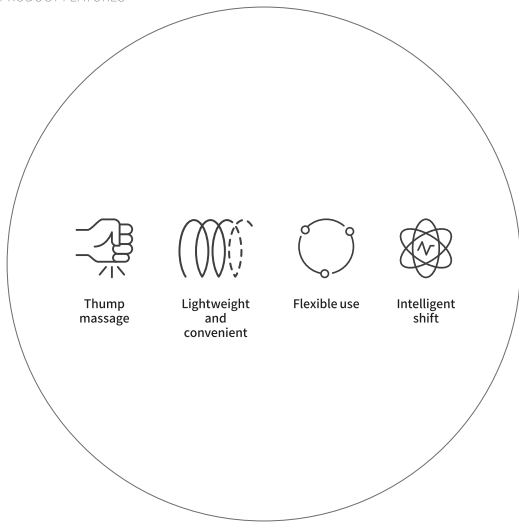


13

## COMPONENT NAMES

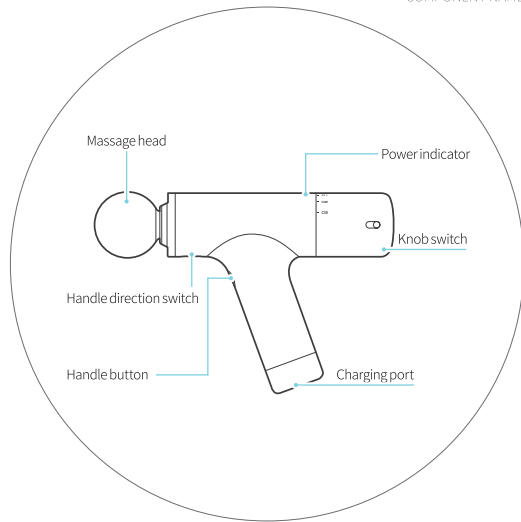

14

## PRODUCT OPERATION

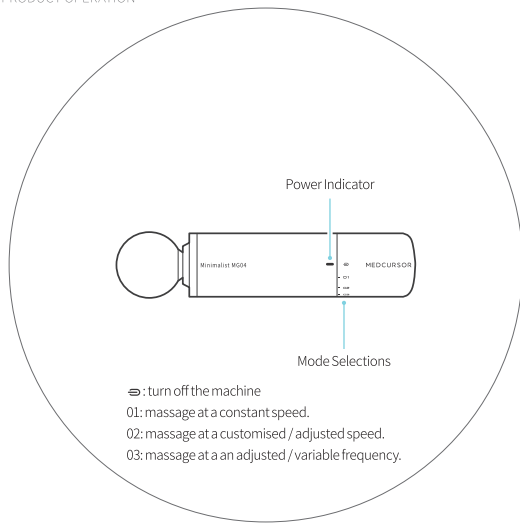

15

EN

## Operation Instructions

### LED Light Indicator

When battery is low, the LED light flashes slowly purple.

The unit automatically stops and the LED light flashes rapidly purple after working for 10 minutes.

When charging, the LED light flashes red.

After being fully charged, the LED light glows green.

When overcurrent/over-temperature protection is functioning, the LED light glows red;

wait until it glows purple, then you may continue to use.

16

EN

### Charging Method

1. Open the charging port cover at the bottom of the handle and connect the charging cable.
2. Plug the adapter into the power socket.

※When charging, the main machine and the adapter may go slightly hot, which is normal.

### Charging Time

1. It takes 2-3 hours to fully charge the device..
2. When fully charged, it supports up to 13 times of use for 10 minutes each. Multiple use in a row is not recommended. Wait 20-30 minutes before next use. (The time interval may vary due to ambient temperature.)
3. While charging, the LED light on the main machine flashes red.
4. Once fully charged, the LED light turns solid green..


### Battery

When battery is low, the LED light starts to flash purple. Please charge before next use.

17

### After Being Fully Charged

Unplug the adapter and disconnect the power cord from the main machine.

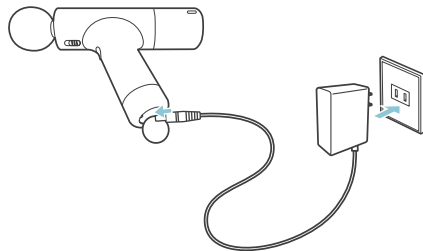

18


EN

### Handle Direction Switch

Push and hold the handle reverse switch located under the front part of the unit (Just above the handle) Do this and then turning or twisting the handle backward/forward (In reverse). (See the image for your reference)

**Note:** The handle can only be reversed up to 180 degrees

To avoid getting your fingers or other body parts pinched, please pause or turn off the machine before conducting directional change.

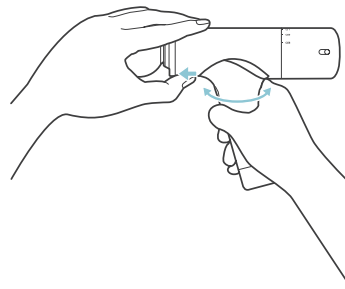

EN

## Proper Usage

Do not use this product while charging. After charging, please unplug the adapter before you use the product.

1. Turn the dial knob switch to turn on the power. The LED light glows and the machine starts working.
2. Turn the dial knob switch to your favorite massage mode or to adjust the massage settings.

### Mode 01

Massage at a constant speed

Under 01, You need to press and hold the handle button for the massager to work.

### Mode 02

Massage at a customised speed

Under 02, the machine will run automatically in this mode. You can change the running speed by holding down the handle button. When it finally adjusted to the speed that you want, release the handle button so the massager will operate at your chosen speed rate.

### Mode 03

Massage at an adjusted and variable frequency

Under 03, You can press and hold the button to increase the shift frequency, the speed changes periodically within a certain range after the button is released. Release the handle button to apply and set the frequency rate that you're in. Note that the speed rate is not fixed.

Double-click the handle button to reset to the initial speed shift so that you can do the adjustment again.

4. The unit automatically stops and the LED light flashes rapidly purple after working for 10minutes. Turn the dial knob switch to OFF position before restarting the machine for continuous use.

5. The product features an overcurrent protection function. For user's safety, the machine stops running automatically and the LED light glows red when excessive vertical load is detected. In this case, move the massage head away and wait for a moment until the LED light turns purple. Then, click the handle button to resume normal use.

6. The handle of the massage gun will slightly heat up when it is working and this is normal. Also, the product has a built-in over heating protection function so there's nothing to worry. Continuous use of this product may lead the motor inside the handle to overheat. Automatic protection will be triggered when the temperature is too high so the unit stops running and the LED light glows red. In this case, stop using it and wait for a moment until LED light turns purple. Then, click the handle button to resume normal use.

#### 7. Turning the Unit off Manually

Turn the dial knob switch to OFF. The unit stops and the LED light goes off. The device shuts down.

#### Note



**Do not disassemble the machine if you are not a professional technician/ electrician.**

## EN

### Troubleshooting

If any faults occur during use, please check the following. If problem still exists or the product doesn't work well, please turn it off and send it to professional and technical staff for repair. Do not disassemble it by yourself!

Basic failures	Possible cause	Solutions
Unable to start	Low battery	Charge the battery in time
Weak force	Low battery	Charge the battery in time
Unit not charging	Adapter not plugged in properly	Check the power adapter if connected to a working power outlet. Connect to another power outlet

If you have any issues, questions or feedback, please send an email to [support@medcursor.com](mailto:support@medcursor.com). We are always happy to help.

## PROTECTION

### EN

#### **Storage**

When the device is stored for more than 3 months, check the device and recharge it before using. Do not store the device closely against the ground, wall or ceiling.

Store in a well-ventilated place and avoid direct sunlight and exposure heat, water or moisture.

#### **Transport**

The product is packed with simple anti-vibration measures for transportation by air, railway, highway and sea. Rain, snow and impact should be avoided. Do not store the unit upside the down.

#### **Power Check**

Check the power supply. Do not start the machine if power supply exceeds the allowed input.

#### **Maintenance**

Keep the unit away from heat, moisture, rodents, dirt and dust. Turn the machine off before cleaning. Store the machine in a safe place after use. Do not place heavy items on the machine. If any fault occurs, call professional and technical staff for repair.

Do not disassemble the machine if you are not a professional technician.

#### **Disposal and Recycling**

This product should not be regarded as a household waste. Please follow local environment laws and guidelines when you have decided to discard this product.

Please remember to respect the local regulations: hand in the non-working electrical equipments to an appropriate waste disposal centre. The packaging material is recyclable. Dispose the packaging in an environmentally friendly manner and make it available for the recyclable material collection-service.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The statements should be displayed in the user manual: changes or modifications not expressly approved by the party

responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:


- reorient or relocate the receiving antenna;
- increase the separation between the equipment and receiver;
- connect the equipment into an outlet on a circuit different from that to which the receiver is connected;
- consult the dealer or an experienced radio/tv technician for

help.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.



EN

Product name	Medcursor Massager Gun
Model	MMG0401
Rated voltage	11.4 V 
Rated power	17 W
Battery capacity	1900mAh
Product weight	615 g
Product size	L170×W150×H48 (mm)

support@medcursor.com

