


# User Manual

Blood Pressure Monitor TMB-1583-BS

Arm Type

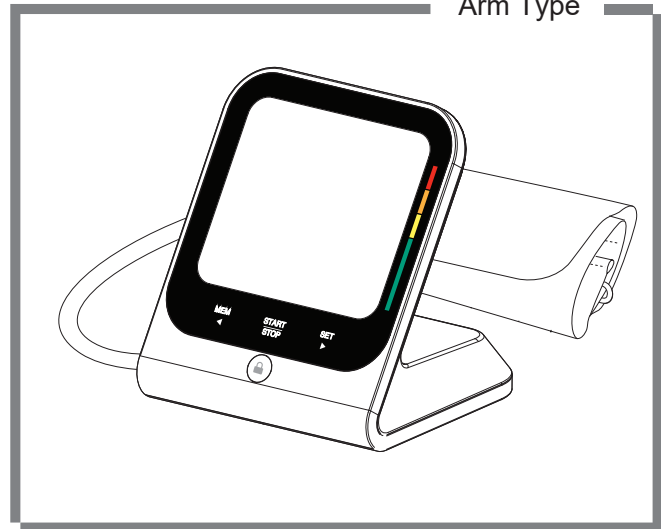

- Thank you for selecting a TRANSTEK Blood Pressure Monitor.
- Please read the user manual carefully and thoroughly to ensure the safe usage of this product. Keep the manual well for further reference in case you have problems.

## Table of Contents

INTRODUCTION.....	2
• General Description	
• Indications for Use	
• Contraindications	
• Measurement Principle	
• Safety Information	
• LCD Display Signal	
• Monitor Components	
• List	
BEFORE YOU START.....	8
• The Choice of Power Supply	
• Installing and Replacing the Batteries	
• Setting Date, Time and Measurement Unit	
• Select the User	
• Pair-up the Blood Pressure Monitor with Your Device	
MEASUREMENT.....	14
• Tie the Cuff	
• Start the Measurement	
DATA MANAGEMENT.....	17
• Recall the Records	
• Delete the Records	
SPECIAL FUNCTION.....	19
• About the Lock Button	
INFORMATION FOR USER.....	20
• Tips for measurement	
• Maintenances	
ABOUT BLOOD PRESSURE.....	22
• What are systolic pressure and diastolic pressure?	
• What is the standard blood pressure classification?	
• Irregular heartbeat detector	
• Why does my blood pressure fluctuate throughout the day?	
• Why do I get a different blood pressure at home compared to the hospital?	
• Is the result the same if measuring on the right arm?	
TROUBLESHOOTING.....	24
SPECIFICATIONS.....	25
AUTHORIZED COMPONENT .....	26
CONTACT INFORMATION.....	26
FCC STATEMENT.....	27
EMC GUIDANCE.....	28

## ♥ General Description

Thank you for selecting TRANSTEK arm type blood pressure monitor (TMB-1583-BS). The monitor features blood pressure measurement, pulse rate measurement and the result storage. The design provides you with two years of reliable service.

Readings taken by the TMB-1583-BS are equivalent to those obtained by a trained observer using the cuff and stethoscope auscultation method.

This manual contains important safety and care information, and provides step by step instructions for using the product.

**Read the manual thoroughly before using the product.**

Features:

- 73 mm×84 mm Digital LCD display
- Maximum 60 records per each user
- 3rd technology: Measuring during inflation

## ♥ Indications for Use

The Transtek Blood Pressure Monitor is a digital monitor intended for use in measuring blood pressure and heartbeat rate with arm circumference ranging from 22 cm ( about 8¾"-12½" ) or 22cm to 42cm(about 8¾"-16½"). It is intended for adult indoor use only.

## ♥ Contraindications






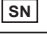





1. The device is not suitable for use on may be pregnant women or pregnant women.
2. The device is not suitable for use on patients with implanted, electrical devices, such as cardiac pacemakers, defibrillators.

## ♥ Measurement Principle

This product uses the Oscillometric Measuring method to detect blood pressure. Before every measurement, the unit establishes a "zero pressure" equivalent to the atmospheric pressure. Then it starts inflating the arm cuff, meanwhile, the unit detects pressure oscillations generated by beat-to-beat pulsatile, which is used to determine the systolic and diastolic pressure, and also pulse rate.

## ♥ Safety Information

The signs below might be in the user manual, labeling or other component. They are the requirement of standard and using.

	Refer to instruction manual/booklet To signify that the instruction manual/booklet must be read.
	Symbol for "Type BF applied part"
	Caution Indicates that caution is necessary when operating the device or control close to where the symbol is placed, or that the current situation needs operator awareness or operator action in order to avoid undesirable consequences.
	The symbol indicates that the product should not be discarded as unsorted waste but must be sent to separate collection facilities for recovery and recycling.
	Symbol for "Recycle"
	Symbol for "Serial Number"
	Symbol for "Direct Current"
	Symbol for "Manufacturer"
	Symbol for "Date and Country of manufacture"
	For indoor use only
	Symbol for "Class II Equipment"


 CAUTION

- \* This device is intended for indoor, home use.
- \* This device is not intended for public use.
- \* This device is portable, but it is not intended for use during patient transport.
- \* This device is not suitable for continuous monitoring during medical emergencies or operations.
- \* This device is intended for non-invasive measuring and monitoring of arterial blood pressure. It is not intended for use on extremities other than the arm, or for any purpose other than obtaining a blood pressure measurement.
- \* This device is for adults. Do not use this device on neonates or infants. Do not use it on children unless otherwise instructed by a medical professional.
- \* Do not use on the women in pregnant, including pre-eclamptic, patients.
- \* The device is not suitable for use on patients with implanted, electrical devices, such as cardiac pacemakers, defibrillators.
- \* The effectiveness of this device has not been established for use:
  - on users with common arrhythmias such as atrial or ventricular premature beats or atrial fibrillation,
  - on users with peripheral arterial disease,
  - on users undergoing intravascular therapy, or with arteriovenous (AV) shunt.
 Consult a medical professional before use.
- \* Do not use this device for diagnosis or treatment of any health problem or disease. Contact your physician if you have or suspect any medical problem. Do not change your medications without the advice of your physician or health care professional.
- \* If you are taking medication, consult your physician to determine the proper time to measure your blood pressure.
- \* This device may be used only for the intended use described in this manual, the manufacturer shall have no liability for any incidental, consequential, or special damages caused by misuse or abuse.
- \* Report any unexpected operation or events to the manufacturer.
- \* Do not apply the cuff on an arm that has an intravenous drip or a blood transfusion attached.
- \* Warning: Do not kink, fold, stretch, compress, or otherwise deform the tube during measuring, as the cuff pressure might continuously increase, which could prevent blood flow and result injury.
- \* Warning: Taking blood pressure measurements too frequently could disrupt blood circulation and cause injuries.
- \* Warning: Do not apply cuff to areas on patient where skin is delicate or damaged. Check cuff site frequently for irritation.
- \* Warning: Do not place the cuff on the arm of a person whose arteries or veins are undergoing medical treatment, i.e. intra-vascular access or intra-vascular therapy or an arteriovenous (A-V) shunt, which could disrupt blood circulation and cause injuries.
- \* Do not place the cuff on the arm on the same side of a mastectomy (especially when lymph nodes have been removed), it is recommended to take measurements on the unaffected side.
- \* Do not wrap the cuff on the same arm to which another monitoring device is applied. One or both devices could temporarily stop functioning if you try to use them at the same time.
- \* Please check that the operation of the device do not result in prolonged impairment of patient blood circulation.
- \* Warning: On the rare occasion of a fault causing the cuff to remain fully inflated during measurement, loosen and remove the cuff immediately. Prolonged high pressure applied to the arm (cuff pressure >300 mmHg or constant pressure >15 mmHg for more than 3 minutes) might lead to bruising and discolored skin.
- \* Warning: Do not use this device with high-frequency (HF) surgical equipment at the same time.

 CAUTION


- \* Warning: This device is not AP/APG equipment. Do not use the device where flammable anesthetic are present, or in environments mixture with air of with oxygen or nitrous oxide.
- \* The device contains sensitive electronic components. To avoid measurement errors, avoid taking blood pressure measurements near a strong electromagnetic field radiated interference signal or electrical fast transient/burst signal.
- \* Wireless communication equipment, such as wireless home network devices, mobile phones, cordless telephones and their base stations, walkie-talkies may cause interference that may affect the accuracy of measurements. A minimum distance of 1 foot (30 cm) should be kept from such devices during a measurement.
- \* You can use this device to take your own measurement, no third-party operator is required.
- \* Please use the device under the environment which is provided in the user manual. Otherwise, the performance and lifetime of the device will be impacted and reduced.
- \* The device may require up to 30 minutes to warm up / cool down from the minimum / maximum storage temperature before it is ready for use.
- \* Warning: Excessive cuff tube lengths could cause strangulation if you don't manage them properly.
- \* Warning: Do not touch output of the batteries/adaptor and the user simultaneously.
- \* Adapter is specified as a part of ME EQUIPMENT.
- \* Warning: The power cord is considered the disconnect device for isolating this equipment from supply mains. Do not position the equipment so that it is difficult to reach or disconnect.
- \* The blood pressure monitor, its adapter, and the cuff are suitable for use within the patient environment.
- \* Warning: Do not use this device if you are allergic to polyester, nylon, or plastic.
- \* Warning: Only use accessories approved by manufacturer. Using unapproved accessories might cause damage to the unit and injure users.
- \* Warning: If you experience discomfort during a measurement, such as pain in the arm or other complaints, press the Power button immediately to release the air from the cuff.
- \* No calibration is required within two years of reliable service.
- \* Do not attempt to repair the unit yourself if it malfunctions. Only have repairs carried out by authorized service centers.
- \* At the request of authorized service personnel, circuit diagrams, component part lists, descriptions, and calibration procedures will be made available by the manufacturer or distributor.
- \* It is recommended that the performance should be checked after repair, maintenance, and every two years of use, by retesting the requirements in limits of the error of the cuff pressure indication and air leakage (testing at least at 50 mmHg and 200 mmHg).
- \* Warning: Do not use the device while under maintenance, or being serviced.
- \* Store your device, cuff and adapter in a clean and dry place, protect it against extreme moisture, heat, lint, dust and direct sunlight. Never place any heavy objects on it.
- \* Make sure the rubber tube of the cuff is not squeezed, stretched, or kinked during storage.
- \* Warning: Keep the device, cuff, and batteries away from children as they may pose a risk of choking or strangulation if used improperly.
- \* Clean both device and cuff with a soft, dry cloth. If necessary use a dampened cloth and natural detergent. Do not use alcohol, benzene, or other harsh chemicals.
- \* Do not wash the cuff in a washing machine or dishwasher!
- \* The service life of the cuff may vary by the frequency of washing, skin condition, and storage state. The typical service life is 10000 times.
- \* Dispose of accessories, detachable parts, and the device according to the local guidelines.

## ♥ LCD display signal

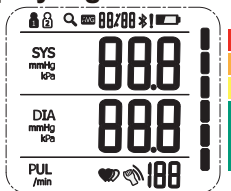


SYMBOL	DESCRIPTION	EXPLANATION
<b>SYS</b>	Systolic blood pressure	High blood pressure
<b>DIA</b>	Diastolic blood pressure	Low blood pressure
<b>PUL</b> /min	Pulse display	Pulse in beats per minute
<b>kPa</b>	kPa	Measurement Unit of the blood pressure (1kPa=7.5mmHg)
<b>mmHg</b>	mmHg	Measurement Unit of the blood pressure (1mmHg=0.133kPa)
<b>Lo +</b>	Low battery	Batteries are low and need to be replaced
<b>♥</b>	Irregular heartbeat	Blood pressure monitor is detecting an irregular heartbeat during measurement.
<b> </b>	Blood pressure level indicator	Indicate the blood pressure level
<b>00:00</b>	Current Time	Month/Day/Year, Hour/Minute
<b>1</b>	User 1	Start measurement,save and transmit the measuring results for User 1
<b>2</b>	User 2	Start measurement,save and transmit the measuring results for User 2
<b>Bluetooth icon</b>	Bluetooth icon	The bluetooth icon blinks when the Bluetooth is working
<b>AVG</b>	The average value	The average value of the latest three records
<b>Motion icon</b>	Motion indicator	Motion may result in an inaccurate measurement
<b>Q</b>	Memory	Indicate it is in the memory mode and which group of memory it is.
<b>!</b>	Bluetooth connection	It blinks when the bluetooth connection fails or the data is not sent.
<b>♥</b>	Heartbeat	Blood pressure monitor is detecting a heartbeat during measurement.

## ♥ Monitor Components

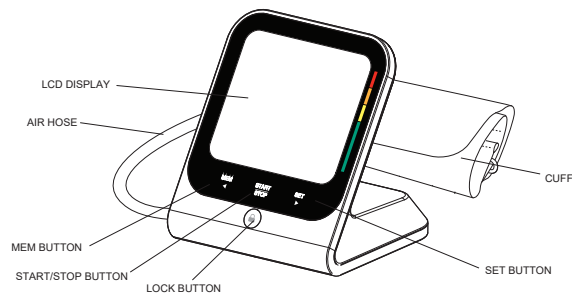

## Component list of pressure measuring system

- 1 Cuff
- 2 Air pipe
- 3 PCBA
- 4 Pump
- 5 Valve

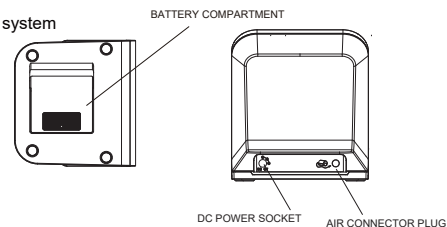

## ♥ List

1. Blood Pressure Monitor (TMB-1583-BS)



3. 4×AAA batteries



2. Cuff (Type BF applied part) (22cm~32cm or 22cm~42cm)




4. User manual

5. AC Adaptor (BLJ06L060100P-U)

## ♥ The Choice of Power Supply

1. Battery powered mode:  
6VDC 4×AAA batteries
2. AC adaptor powered mode:  
6V  $\overline{\rightleftharpoons}$  1A  
(Please only use the recommended AC adaptor model).

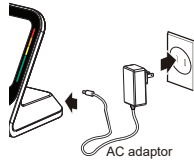


Please unplug the adaptor to depart from the using utility power.

### ⚠ CAUTION

In order to get the best effect and protect your monitor, please use the right battery and special power adaptor which complies with local safety standard.

## ♥ Installing and Replacing the Batteries

- Open the battery cover.
- Install the batteries by matching the correct polarity, as shown.
- Close the cover.

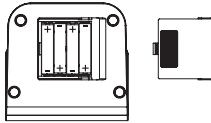

The typical service life of the new and unused batteries is 300 measurements for the operation time is 60s.

Replace the batteries whenever the below happen

- The  $\text{LO} + \square$  shows
- The display is dim.
- The display does not light up

### ⚠ CAUTION

- Do not use new and used batteries together.
- Do not use different types of batteries together.
- Do not dispose the batteries in fire. Batteries may explode or leak.
- Remove batteries if the device is not likely to be used for some time.
- Worn batteries are harmful to the environment. Do not dispose with daily garbage.
- Remove the old batteries from the device following your local recycling guidelines.

## ♥ Setting Date, Time and Measurement Unit

It is important to set the clock before using your blood pressure monitor, so that a time stamp can be assigned to each record that is stored in the memory. (The setting range of the year :2018—2058 time format:24H)

1. When the monitor is off, hold pressing "SET" button for 3 seconds to enter the mode for year setting.




2. Press the "MEM" button to change the [YEAR]. Each press will increase the numeral by one in a cycling manner.




3. When you get the right year, press "SET" button to confirm then it will turn to next step.




4.Repeat step 2 and 3 to set the [MONTH] and [DAY].

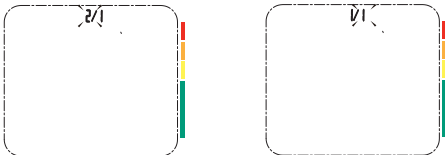

5.Repeat step 2 and 3 to set the [HOUR] and [MINUTE].

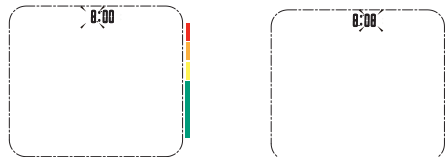

6.Repeat step 2 and 3 to set the [UNIT].

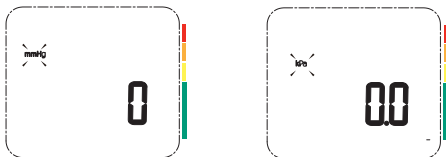

7 After UNIT is set,the LCD will display "donE" first, then display all the settings you have done and then turn off.



## ♥ Select the User


1.When the monitor is off , press and hold the MEM button to enter user setting mode. The user ID will blink.



2.Then press MEM button again, select the user ID between user 1 and user 2.



3. After selecting the suitable user ID, press SET button to confirm. The LCD will show " User ID+ donE " and then turn off.

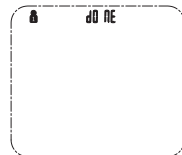

## ▼ Pair-up the Blood Pressure Monitor with Your Device

1. Turn on Bluetooth and the app. Make sure both are ON when pair-up is proceeding.
2. When the monitor is OFF, press and hold the START/STOP button to start pair-up. The symbol "⌘" and user ID will be shown on the LCD alternatively, indicating pair-up is proceeding.



If **SUCCEED**, the symbol "⌘ + [ ]" and user ID will be shown on the LCD.



If **FAIL**, the symbol "⌘ + E 12" and user ID will be shown on the LCD.



3. The monitor will shut off after Pair-up process is complete.

**Bluetooth Module No.:** LS51802  
**RF Frequency Range:** 2402 MHz to 2480 MHz  
**Output Power Range:** -2.6 dBm  
**Supply Voltage:** 2V-3.6 V  
**Transmitting Distance:** 10 meters

List of compatible devices:

For iOS devices:

The operating system must be iOS 8 or more, such as iPhone 4S, iPhone 5/5C/5S, iPhone 6/6 Plus and so on.

For Android devices:

The operating system must be 4.3 or more.

### ⚠ CAUTION

- Interference may occur in the vicinity of equipment marked with the following symbol (⚠). And TMB-1583-BS may interfering vicinity electrical equipment.
- Sensitive people, including pregnant women pre-eclamptic and those who implanted medical electronic instruments, should avoid using the unit whenever possible.
- Keep the monitor at least 20 centimeters away from the human body (especially the head) when the data transmission is proceeding after measurement.
- To enable the data transmission function, this product should be paired to Bluetooth end at 2.4 GHz.


#### How to mitigate possible interference?

1. The range between the device and BT end should be reasonably close, from 1 meter to 10 meters. Please ensure no obstacles between the device and BT end so as to obtain quality connection and to lower the RF output range.
2. To avoid interference, other electronic devices (particularly those with wireless transmission / Transmitter) should be kept at least 1 meter away from the monitor.



## ♥ Tie the cuff

1. Tie the cuff on your upper arm, then position the tube off-center toward the inner side of arm in line with the little finger. Or position the artery mark  $\Phi$  over the main artery (on the inside of your arm). Note: Locate the main artery by pressing with 2 fingers approximately 2 cm above the bend of your elbow on the inside of your left arm. Identify where the pulse can be felt the strongest. This is your main artery.

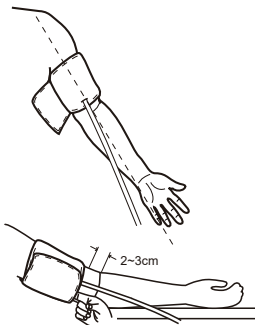


2. The cuff should be snug but not too tight. You should be able to insert one finger between the cuff and your arm.

3. Sit comfortably with your tested arm resting on a flat surface.

### 4. Patients with Hypertension:

The middle of the cuff should be at the level of the right atrium of the heart; Before starting measurement, please sit comfortably with legs uncrossed, feet flat on the floor, back and arm supported.


- Rest for 5 minutes before measuring.
- Wait at least 3 minutes between measurements. This allows your blood circulation to recover.
- For a meaningful comparison, try to measure under similar conditions. For example, take daily measurements at approximately the same time, position of upper arm, or as directed by a physician.

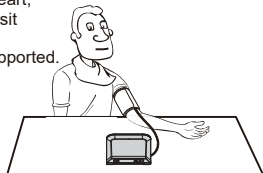


## ♥ Start the Measurement

1. When the monitor is off, press the "START/STOP" button to turn on the monitor, and it will finish the whole measurement. (Take user 1 for example)

LCD display

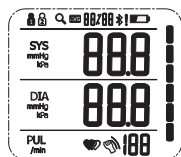


Adjust the zero.

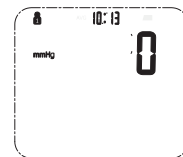


Inflating and measuring.




Display and save the measurement results.

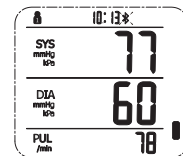


2. This device will proceed to data transmission after measurement. The Bluetooth symbol blinks on the LCD indicates data is transmitting.

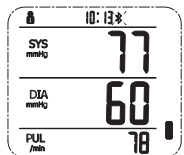

If the data transmission fails, the LCD will display "1".

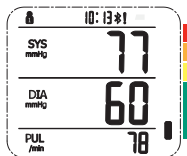


If the data transmission succeeds, the symbol "✕" will disappear.

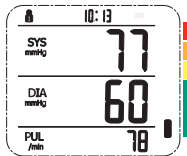


3. Press the "START/STOP" button to power off, otherwise it will turn off within 1 minute.

## ♥ Recall the Records

1. When the monitor is off, please press the "MEM" to show the average value of the latest three records. If the records are less than three groups, it will display the latest record first. (Take user 1 for example)

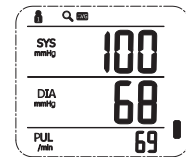

2. Press "MEM" button or "SET" button to get the record you want.

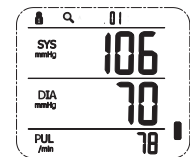

The current No., date and time of the record will be shown alternately.

### ⚠ CAUTION

The most recent record (1) is shown first. Each new measurement is assigned to the first (1) record. All other records are pushed back one digit (e.g., 2 becomes 3, and so on), and the last record (60) is dropped from the list.

## ♥ Delete the Records

If you did not get the correct measurement, you can delete all results by following steps below. (Take user 1 for example)

1. Hold pressing "MEM" button for 3 seconds when the monitor is in the memory recall mode, the "User ID + dEL ALL" will flash on the display.




2. Press "SET" button to confirm deleting and the monitor will display "dEL dOnE"+User ID and then turn off.



Note: To exit the delete mode without deleting any records, press START/STOP button before pressing "SET" to confirm any delete commands.

3. If there is no record, the following display will be shown.

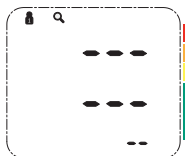

## ♥ About the Lock Button

The unnecessary keys' touch will make the blood pressure monitor turn on and waste electricity. To avoid this, you can press the Lock button to lock the keys if necessary.

Hold pressing the Lock button until the LCD displays "OFF", this indicates the touch keys (such as MEM, START/STOP, SET) have been locked.



Hold pressing the Lock button until the LCD displays "ON" to unlock the keys. You can press any buttons to use the blood pressure monitor.



## ♥ Tips for Measurement

Measurements may be inaccurate if taken in the following circumstances.



Within 1 hour  
after dinner or drinking



Immediate measurement  
after tea, coffee, smoking



Within 20 minutes  
after taking a bath



When talking or moving your fingers



In a very cold environment



When you want to discharge urine




## ♥ Maintenance

In order to get the best performance, please follow the instructions below.



Put in a dry place and avoid the sunshine

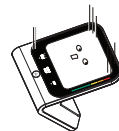

Avoid intense shaking  
and collisions



Using wet cloths to remove dirt



Avoid touching water,  
clean it with a dry cloth in case.




Avoid dusty and unstable  
temperature environment



Do not attempt to clean the reusable cuff  
with water and never immerse the cuff in  
water.

### ♥ What are systolic pressure and diastolic pressure?

When ventricles contract and pump blood out of the heart, the blood pressure reaches its maximum value in the cycle, which is called systolic pressure. When the ventricles relax, the blood pressure reaches its minimum value in the cycle, which is called diastolic pressure.

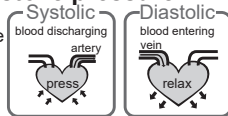

### ♥ What is the standard blood pressure classification?

This chart reflects blood pressure categories defined by American Heart Association.

Blood Pressure Category	Systolic mmHg (upper#)		Diastolic mmHg (lower#)
Normal	less than 120	and	less than 80
Elevated	120-129	and	less than 80
High Blood Pressure (Hypertension) Stage 1	130-139	or	80-89
High Blood Pressure (Hypertension) Stage 2	140 or higher	or	90 or higher
Hypertensive Crisis (Consult your doctor immediately)	Higher than 180	and/or	Higher than 120

The chart on the right is the standard blood pressure classification published by American Heart Association (AHA).

#### ⚠ CAUTION

Please consult a physician if your measuring result falls outside the range. Please note that only a physician can tell whether your blood pressure value has reached a dangerous point.

### ♥ Irregular Heartbeat Detector

An irregular heartbeat is detected when a heartbeat rhythm varies while the device is measuring systolic pressure and diastolic pressure. During each measurement, blood pressure monitor will keep a record of all the pulse intervals and calculate the average value of them. If there are two or more pulse intervals, the difference between each interval and the average is more than the average value of  $\pm 25\%$ , or there are four or more pulse intervals, the difference between each interval and the average is more than the average value of  $\pm 15\%$ , then the irregular heartbeat symbol will appear on the display with the measurement result.

#### ⚠ CAUTION

The appearance of the IHB icon indicates that a pulse irregularity consistent with an irregular heart-beat was detected during measurement. Usually this is NOT a cause for concern. However, if the symbol appears often, we recommend you seek medical advice. Please note that the device does not replace a cardiac examination, but serves to detect pulse irregularities at an early stage.

### ♥ Why does my blood pressure fluctuate throughout the day?

1. Individual blood pressure varies many times everyday. It is also affected by the way you tie your cuff and your measurement position, so please take the measurement under the same conditions.
2. If the person takes medicine, the pressure will vary more.
3. Wait at least 3 minutes for another measurement.



### ♥ Why do I get a different blood pressure at home compared to the hospital?

The blood pressure is different even throughout the day due to weather, emotion, exercise etc. Also, there is the "white coat" effect, which means blood pressure usually increases in clinical settings.


What you need to pay attention to when you measure your blood pressure at home: If the cuff is tied properly. If the cuff is too tight or too loose. If the cuff is tied on the upper arm. If you feel anxious. Taking 2-3 deep breaths before beginning will be better for measuring. Advice: Relax yourself for 4-5 minutes until you calm down.


### ♥ Is the result the same if measuring on the right arm?

It is ok for both arms, but there will be some different results for different people. We suggest you measure the same arm every time.



This section includes a list of error messages and frequently asked questions for problems you may encounter with your blood pressure monitor. If the products not operating as you think it should, check here before arranging for servicing.


PROBLEM	SYMPTOM	CHECK THIS	REMEDY
No power	Display will not light up.	Batteries are exhausted.	Replace with new batteries
		Batteries are inserted incorrectly.	Insert the batteries correctly
		AC adaptor is inserted incorrectly.	Insert the AC adaptor tightly
Low batteries	Display is dim or show 	Batteries are low.	Replace with new batteries
Error message	E 01 shows	The cuff is too tight or too loose.	Refasten the cuff and then measure again.
	E 02 shows	The monitor detected motion, talking, or the pulse is too poor while measuring.	Relax for a moment and then measure again.
	E 03 shows	The measurement process does not detect the pulse signal.	Loosen the clothing on the arm and then measure again.
	E 04 shows	The treatment of the measurement failed.	Relax for a moment and then measure again.
	EExx, shows on the display.	A calibration error occurred.	Retake the measurement. If the problem persists, contact the retailer or our customer service department for further assistance. Refer to the warranty for contact information and return instructions.
Warning message	"out " shows	Out of measurement range	Relax for a moment. Refasten the cuff and then measure again. If the problem persists, contact your physician.

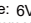Power supply	Battery powered mode: 6VDC 4×AAA batteries AC adaptor powered mode: 6V  1A (Please only use the recommended AC adaptor model).
Display mode	Digital LCD V.A.84mm×73mm
Measurement mode	Oscillographic testing mode
Measurement range	Rated cuff pressure: 0mmHg~299mmHg(0kPa ~ 39.9kPa) Measurement pressure: SYS: 60mmHg~230mmHg (8.0kPa~30.7kPa) DIA: 40mmHg~130mmHg (5.3kPa~17.3kPa) Pulse value: (40-199)beat/minute
Accuracy	Pressure: 5°C~40°C with $\pm 0.4$ kPa(3mmHg) Pulse value: $\pm 5\%$
Normal working condition	A temperature range of :+5°C to +40°C A relative humidity range of 15% to 90%, non-condensing, but not requiring a water vapour partial pressure greater than 50 hPa An atmospheric pressure range of : 700 hPa to 1060 hPa
Storage & transportation condition	Temperature:-20°C to +60°C A relative humidity range of $\leq 93\%$ , non-condensing, at a water vapour pressure up to 50hPa
Measurement perimeter of the upper arm	About 22cm~32cm or 22cm~42cm
Net Weight	Approx.250g(Excluding the dry cells and cuff)
External dimensions	Approx.107mm×103mm×118mm
Attachment	4×AAA batteries, user manual , AC adaptor
Mode of operation	Continuous operation
Degree of protection	Type BF applied part
Protection against ingress of water	IP21 It means the device could protected against solid foreign objects of 12.5mm and greater, and protect against vertically falling water drops.
Device Classification	Battery Powered Mode: Internally Powered ME Equipment AC Adaptor Powered Mode: Class II ME Equipment
Software Version	A01

WARNING: No modification of this equipment is allowed.

## ♥ Authorized Component

1. please use the TRANSTEK authorized adapter.

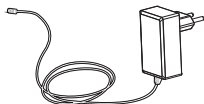

Adapter  
 Model:BLJ06L060100P-U  
 Input: AC 100-240V  
 50/60Hz 0.2A Max  
 Output: 6V $\overline{\text{---}}$  1000mA

## ♥ Contact Information

For more information about our products, please visit [www.transtekcorp.com](http://www.transtekcorp.com). you can get customer service, usual problems and customer download, transtek will serve you anytime.

**Manufactured by:** Guangdong Transtek Medical Electronics Co., Ltd.  
**Company:** Guangdong Transtek Medical Electronics Co., Ltd.  
**Address:** Zone A, No.105, Dongli Road, Torch Development District,  
 528437 Zhongshan, Guangdong, China

## ♥ FCC Statement

FCC ID:OU9TMB1583B2

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Caution: The user is cautioned that changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC Radiation Exposure Statement:

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

## ♥ EMC Guidance

The ME EQUIPMENT or ME SYSTEM is suitable for home healthcare environments.

Warning: Don't be near the active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.

Warning: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

Warning: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

Warning: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the equipment TMB-1583-BS including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

Technical description:

- All necessary instructions for maintaining BASIC SAFETY and ESSENTIAL PERFORMANCE with regard to electromagnetic disturbances for the excepted service life.
- Guidance and manufacturer's declaration-electromagnetic emissions and Immunity.

Table 1

Guidance and manufacturer's declaration - electromagnetic emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class [ B ]
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations / flicker emissions IEC 61000-3-3	Comply

Table 2

Guidance and manufacturer's declaration – electromagnetic Immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV signal input/output 100 kHz repetition frequency	±2 kV for power supply lines Not Applicable 100 kHz repetition frequency
Surge IEC61000-4-5	±0.5 kV, ±1 kV differential mode ±0.5 kV, ±1 kV, ±2 kV common mode	±0.5 kV, ±1 kV differential mode Not Applicable
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% U <sub>T</sub> ; 0.5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°. 0% U <sub>T</sub> ; 1 cycle and 70% U <sub>T</sub> ; 25/30 cycles; Single phase: at 0°. 0% U <sub>T</sub> ; 250 / 300 cycle	0% U <sub>T</sub> ; 0.5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°. 0% U <sub>T</sub> ; 1 cycle and 70% U <sub>T</sub> ; 25/30 cycles; Single phase: at 0°. 0% U <sub>T</sub> ; 250 / 300 cycle
Power frequency magnetic field IEC 61000-4-8	30 A/m 50 Hz / 60 Hz	30 A/m 50 Hz / 60 Hz
Conducted RF IEC61000-4-6	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80% AM at 1 kHz	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80% AM at 1 kHz
Radiated RF IEC61000-4-3	10 V/m 80 MHz – 2,7 GHz 80% AM at 1 kHz	10 V/m 80 MHz – 2,7 GHz 80% AM at 1 kHz
NOTE U <sub>T</sub> is the a.c. mains voltage prior to application of the test level.		



## ♥ EMC Guidance

The ME EQUIPMENT or ME SYSTEM is suitable for home healthcare environments.

Warning: Don't be near the active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.

Warning: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

Warning: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

Warning: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the equipment TMB-1583-BS including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

Technical description:

- All necessary instructions for maintaining BASIC SAFETY and ESSENTIAL PERFORMANCE with regard to electromagnetic disturbances for the excepted service life.
- Guidance and manufacturer's declaration-electromagnetic emissions and Immunity.

Table 1

Guidance and manufacturer's declaration - electromagnetic emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class [ B ]
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations / flicker emissions IEC 61000-3-3	Comply

Table 2

Guidance and manufacturer's declaration – electromagnetic Immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV signal input/output 100 kHz repetition frequency	±2 kV for power supply lines Not Applicable 100 kHz repetition frequency
Surge IEC61000-4-5	±0.5 kV, ±1 kV differential mode ±0.5 kV, ±1 kV, ±2 kV common mode	±0.5 kV, ±1 kV differential mode Not Applicable
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% U <sub>T</sub> ; 0.5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°. 0% U <sub>T</sub> ; 1 cycle and 70% U <sub>T</sub> ; 25/30 cycles; Single phase: at 0°. 0% U <sub>T</sub> ; 250 / 300 cycle	0% U <sub>T</sub> ; 0.5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°. 0% U <sub>T</sub> ; 1 cycle and 70% U <sub>T</sub> ; 25/30 cycles; Single phase: at 0°. 0% U <sub>T</sub> ; 250 / 300 cycle
Power frequency magnetic field IEC 61000-4-8	30 A/m 50 Hz / 60 Hz	30 A/m 50 Hz / 60 Hz
Conducted RF IEC61000-4-6	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80% AM at 1 kHz	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80% AM at 1 kHz
Radiated RF IEC61000-4-3	10 V/m 80 MHz – 2,7 GHz 80% AM at 1 kHz	10 V/m 80 MHz – 2,7 GHz 80% AM at 1 kHz
NOTE U <sub>T</sub> is the a.c. mains voltage prior to application of the test level.		

Table 3

Guidance and manufacturer's declaration - electromagnetic immunity								
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment)	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Maximum Power (W)	Distance (m)	IEC 60601-1-2 Test Level (V/m)	Compliance level (V/m)
	385	380-390	TETRA 400	Pulse modulation 18 Hz	1.8	0.3	27	27
	450	430-470	GMRS 460, FRS 460	FM $\pm$ 5k Hz deviation 1 kHz sine	2	0.3	28	28
	710	704-787	LTE Band 13, 17	Pulse modulation 217 Hz	0.2	0.3	9	9
	745							
	780							
	810	800-960	GSM 800/900, TETRA 800, IDEN 820, CDMA 850, LTE Band 5	Pulse modulation 18 Hz	2	0.3	28	28
	870							
	930							
	1720							
	1845	1700-1990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation 217 Hz	2	0.3	28	28
	1970							
	2450	2400-2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0.3	28	28
	5240	5100-5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	0.2	0.3	9	9
	5500							
	5785							
	5785							